


making a difference

PREANALYTICS CATALOG



making a difference

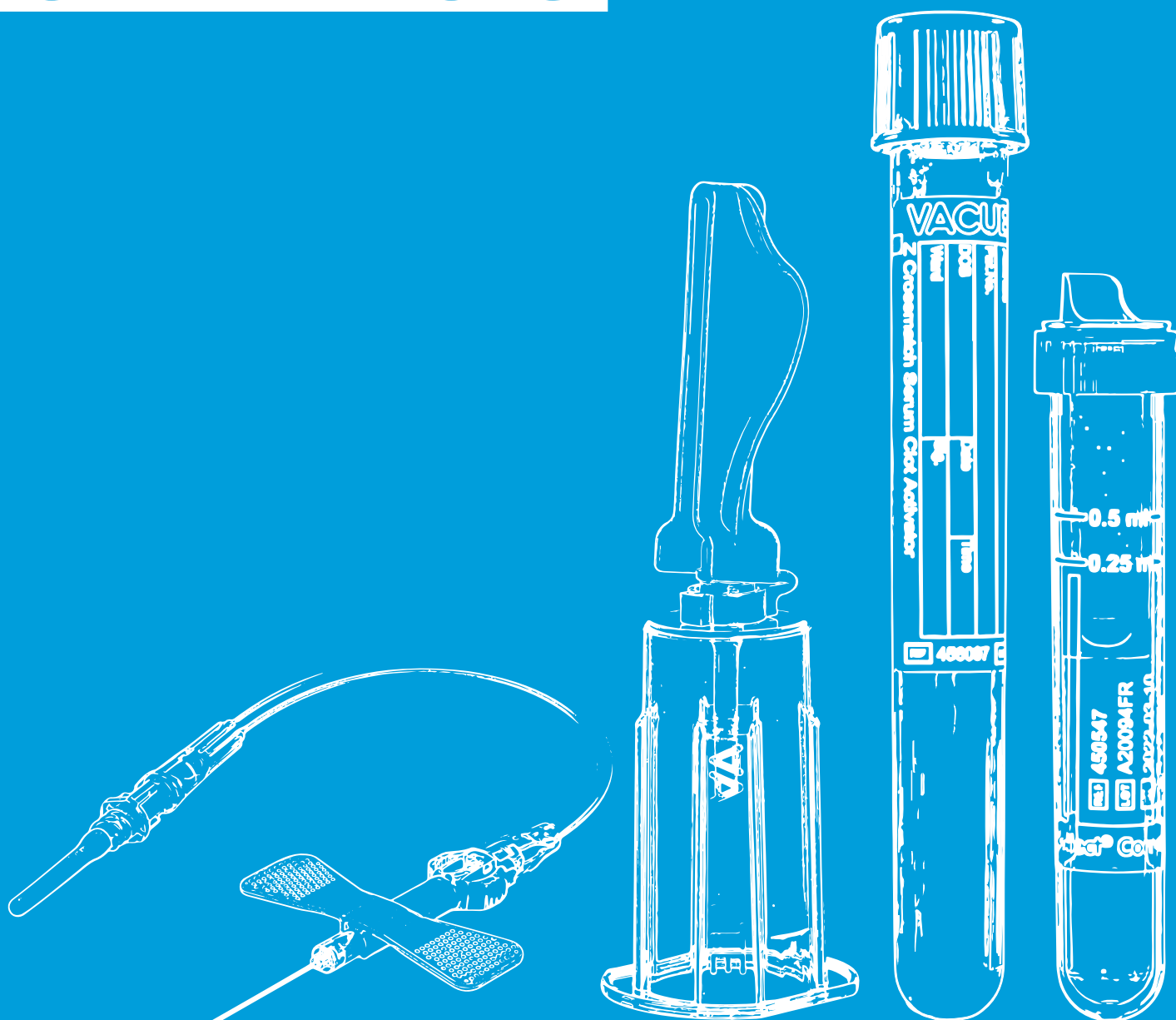

To us, health is the highest personal priority.
This inspires our daily mission of serving life sciences
and healthcare professionals with advanced and
sustainable products and solutions.

It is the passion and commitment of our employees
employees that makes the difference.
With a targeted approach, our global team
continuously strives for excellence
and works closely with our customers.

We empower people to make a difference in health.

CONTENT


6



PRODUCTS

Here you can find
our product
overview


10

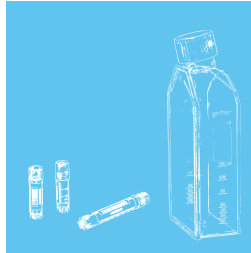

BIOSCIENCE

Making a differ-
ence in research
and science


12

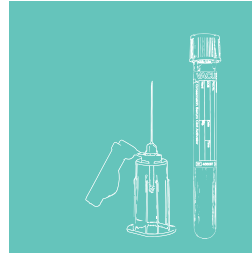

PREANALYTICS

Prenalytical
specimen collec-
tion systems and
safety solutions


14

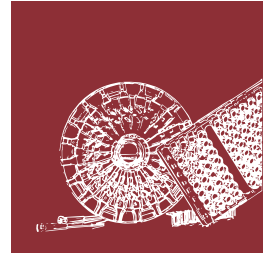

OEM

Tailor made
solutions

17

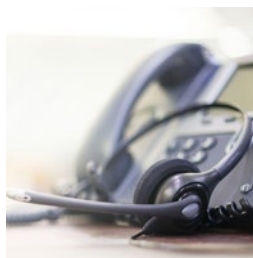

SERVICE

E-Commerce
Solutions from
Greiner Bio-One

18



VALUE PROPOSITIONS

About our value
propositions

284



NUMERICAL INDEX

35

VACUETTE® BLOOD COLLECTION TUBES

VACUETTE® Coagulation Tubes	40
VACUETTE® Serum Tubes	43
VACUETTE® Heparin Tubes	48
VACUETTE® EDTA Tubes	52
VACUETTE® No Additive Tubes	56
VACUETTE® Glycolytic Inhibitor Tubes	58
VACUETTE® Trace Element Tubes	61
VACUETTE® Crossmatch Tubes	63

67

MINICOLLECT® BLOOD COLLECTION TUBES

MiniCollect® Coagulation Tubes	72
MiniCollect® Serum Tubes	74
MiniCollect® Heparin Tubes	76
MiniCollect® EDTA Tubes	78
MiniCollect® Glycolytic Inhibitor Tubes	80

83

SAFETY PRODUCTS

VACUETTE® QUICKSHIELD	88
VACUETTE® EVOPROTECT	91
SAFETY Blood Collection Set	94
MiniCollect® Safety Lancets	97

101

IV CATHETERS

CLiP® Automatic Safety IV Catheter	102
------------------------------------	-----



shop.gbo.com

Online Shop

107

SHARPS

Blood Collection Needles	108
Blood Collection Sets	110

113

BLOOD COLLECTION ACCESSORIES

Tube Holder	114
VACUETTE® SECONDARY Tubes	118
Luer Adapter	120
VACUETTE® Tourniquets	122
VACUETTE® ID Snap Rings	125
Blood Collection Tray	127

131

DISPOSAL

Disposal Containers	132
---------------------	-----

137

SAMPLE TRANSPORT

VACUETTE® Transport Line	138
VACUETTE® Transport Line Accessories	140

143

GENERAL LABORATORY EQUIPMENT

Snap Caps	144
MiniCollect® Accessories Tubes	146

149

DEVICES

VeinViewer®	150
-------------	-----

153

URINE COLLECTION SYSTEM

VACUETTE® Urine Tubes	154
Urine Sets	157
Accessories Urine Collection	160

165

MOLECULAR BIOLOGY



VACUETTE® Virus Stabilization Tube	166
------------------------------------	-----



shop.gbo.com

Online Shop

BIO SCIENCE

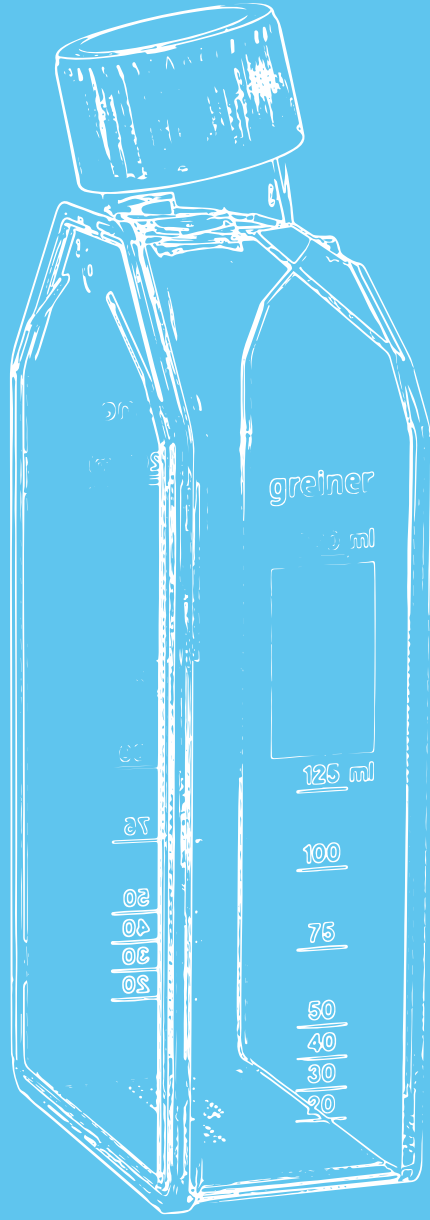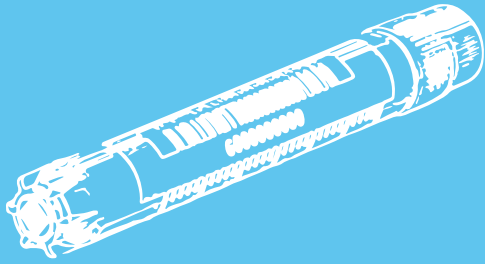

MAKING A DIFFERENCE IN RESEARCH AND SCIENCE

Greiner Bio-One is one of the leading suppliers of special products for the cultivation and analysis of cell and tissue cultures. Under the brand name CELLSTAR®, we offer cell culture flasks, dishes and plates in a wide variety of formats and surface modifications, so that you will always find exactly the right product for every application and all cell types.


In addition, we have been developing and producing microplates for high-throughput screening for more than 60 years, enabling

industry and research to carry out fast and efficient drug testing. Among other things, Greiner Bio-One was the first manufacturer to introduce a microplate in 1536 well format to the market, setting new standards in terms of automation, performance and cost reduction.

Based on decades of experience in the cryogenic storage of samples, we also offer comprehensive solutions for automated storage systems in biobanks.

- / High-quality products for your lab
- / Special products for cell and tissue culture
- / Microplates for industry and research
- / Cryogenic storage systems

PREANALYTICS

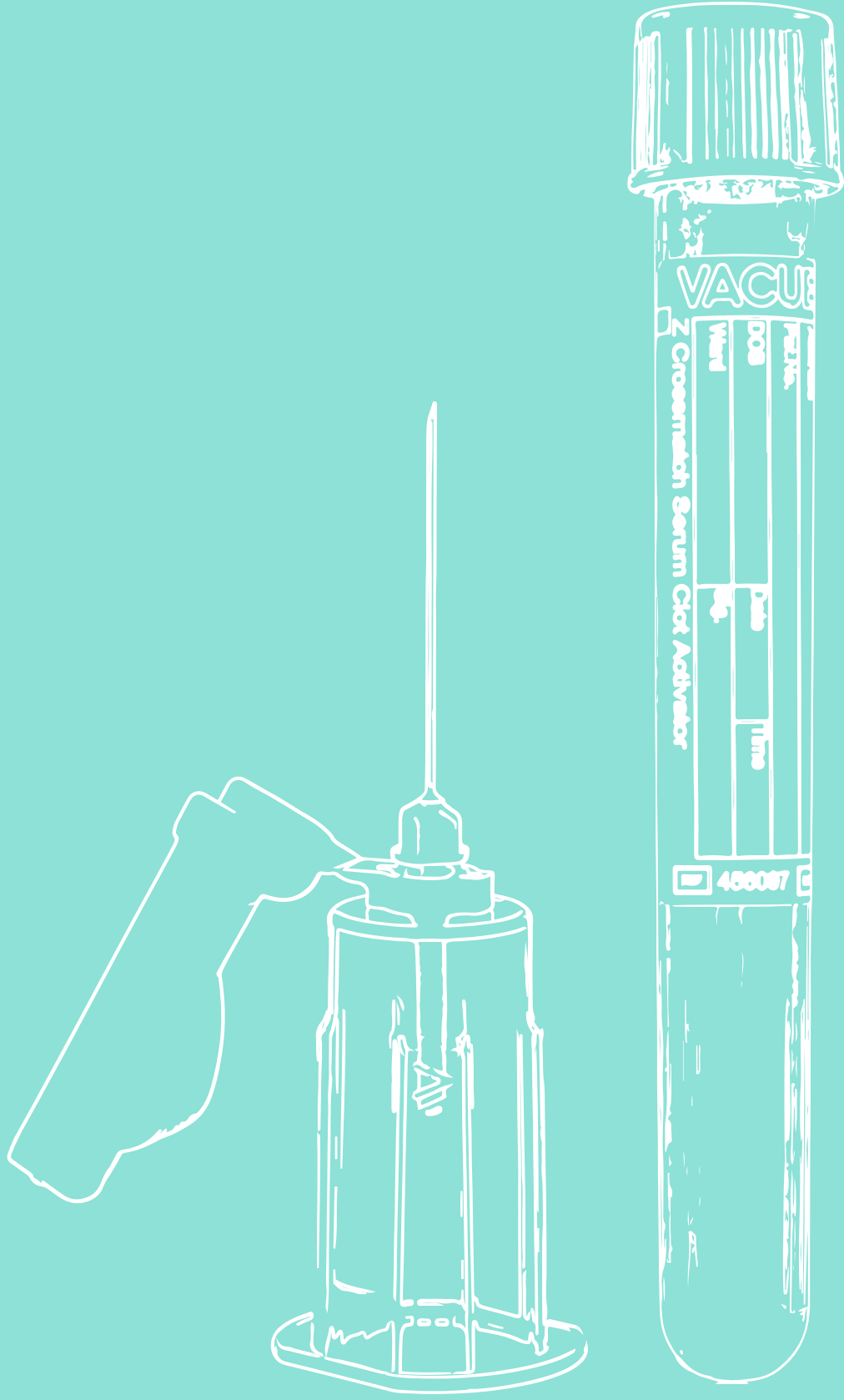

PREANALYTICAL SPECIMEN COLLECTION SYSTEMS AND SAFETY SOLUTIONS

Greiner Bio-One was the first company in the world to successfully launch a blood collection system made out of virtually unbreakable PET plastic including a safety cap. Many components of the VACUETTE® System have often been copied, but never reached our high quality standards.


As a trendsetter in the field

of sample collection, we offer solutions that meet the highest demands of our customers. Because your work is an important contribution to health protection.

Continuously new and improved product solutions, perfect functionality, maximum safety and high product quality are the result of decades of experience.

- / Provide safety
- / Simplify handling
- / Provide support
- / Increase efficiency
- / Save costs

OEM

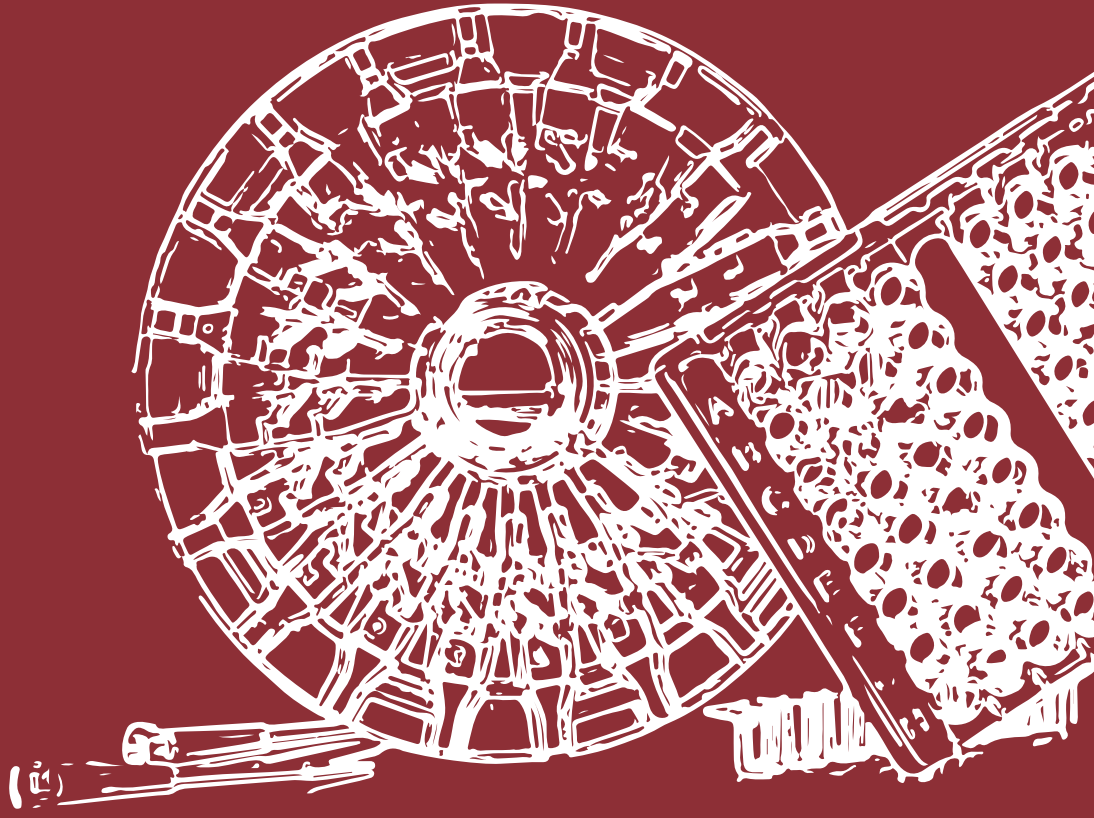

TAILOR MADE SOLUTIONS

Greiner Bio-One as an Original Equipment Manufacturer (OEM) is a long-term partner of the pharmaceutical industry, biotechnology, diagnostic and medical technology industries. The company manufactures numerous products through injection moulding and offers the entire product development and production process. From design,

via prototyping, through fully automated manufacturing, Greiner Bio-One delivers the complete solution with support of experienced scientists, engineers and specialists.

With a global sales network specifically established for the field of OEM, Greiner Bio-One offers custom-made solutions with personal service.

- / Research & development with diverse application and research laboratories**
- / Construction of moulds, planning of plants and assembly lines**
- / Modern, fully automated production**
- / Post-production according to customer requirements**
- / Comprehensive quality assurance and management system**
- / Worldwide networked warehousing and customer service**



FAST. SECURE. EFFICIENT. E-COMMERCE SOLUTIONS FROM GREINER BIO-ONE

As your strong partner in the fields of biotechnology, diagnostics, medical devices and research, Greiner Bio-One already provides you with solutions on a daily basis. In addition, we can offer you a wide selection of e-business solutions to ensure the optimal integration of our products and services into your purchasing process. That in turn will significantly reduce your process-related costs and warehouse costs.

Take advantage of e-business connectivity and make your processes faster, more secure and more efficient.

Greiner Bio-One offers you numerous solutions and the necessary expertise to support your processes, regardless of which ERP system you use. Would you like to purchase via our online shop or are you interested in one of our eProcurement solutions such as EDI or OCI?


In order to stay informed, you can subscribe to our Greiner Bio-One Newsletter and never miss any news about our products and events.



- / Online-Shop
- / EDI (Electronic Data Interchange)
- / OCI (Open Catalog Interface)
- / Electronic catalogs



Subscribe to our newsletter
and never miss
the latest news!

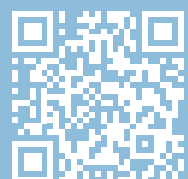



As scientists, we
share your passion for
your research success.

We are happy to support you
wherever and whenever
we are needed.

TRUE TRUST

Take advantage of our experience and passion for your projects and tasks



WE KEEP OUR PROMISES

For us, real trust means keeping promises. This is one of the core values that we are committed to as a company, both towards you as our customers and towards our employees.

As a family-owned company, we believe that trusting, active partnership is the key to sustainable and reliable processes. We

understand your requirements, speak your language and are happy to support you in advancing your scientific progress. We are at your side at every stage – our competent team of experts is always on hand, from the first contact and flawless delivery of our products, to all forms of after-sales service.

PUT YOUR TRUST IN US:

- / Partner for sustainable and reliable processes
- / More than 60 years of experience in science and healthcare
- / Innovation with tradition
- / Support during every work step



DISCOVER
MORE:



With state-of-the-art technology,
we develop high-performance
solutions for consistent,
reproducible results that
you can always rely on.



QUALITY

Reproducible results
you can rely on

CONSISTENTLY HIGH QUALITY FOR RELIABLE RESULTS

Quality has been our top priority for more than 60 years. This applies all over the world, as all of our production facilities are committed to the same high standard.

Our products are made from strictly controlled, certified raw materials. With state-of-the-art production technology and

continuous in-process checks by specially trained employees, we ensure consistent quality of the products. Quality promise not only relates to the manufacturing process and therefore to product quality, but to every process and every interaction with you. We are ISO 9001 and ISO 13485 certified.



CONSISTENT HIGH QUALITY:

- / Location-independent
- / State-of-the-art manufacturing technologies
- / Controlled and certified raw materials
- / ISO 9001 and ISO 13485 certified



DISCOVER MORE:





Thanks to our presence in over 100
countries, you can rely on us to
provide you with consistent high
quality, availability and support.
Wherever you are.

GLOBAL PRESENCE

There to support you
- wherever you are

ON-SITE FOR YOU WORLDWIDE

By global presence we mean being there for you whenever and wherever you need us.

Our seven production sites in Europe, Asia and North and South America meet the same high quality standards and

ensure low-risk delivery routes and a secure product supply. The support teams of our 28 offices and distributors in more than 100 countries are helpful and competent, and are always on hand when you need help or support.



WE PLACE IMPORTANCE ON:

- / **Competent local support**
- / **Optimum product availability**
- / **Consistent high quality worldwide**



**DISCOVER
MORE:**



A man with a shaved head and glasses, wearing a light blue button-down shirt, is standing in front of a glass wall or window. He is holding a black marker in his right hand and has just finished drawing a lightbulb on the glass. The lightbulb is drawn with simple black lines, showing the filament and the base. The background is a blurred office interior with wooden paneling and other office furniture.

Our innovative spirit
supports the results of
your pioneering work.
In the past, present
and future.

IDEAS BECOME SOLUTIONS

For more than 60 years, our passion for plastic technology and manufacturing processes has driven us to constantly innovate to support you optimally in your daily work and ensure your success.

We have been instrumental in the development of consumables in the diagnostics and life sciences industries, as well as in the medical market, since 1963.

We are proud of our rich pioneering history – from the production of the first plastic petri dishes, to Europe's first plastic vacuum sampling system or the invention of the world's first 1536-well microplates. We will continue to strive for new and innovative product solutions that are precisely tailored to your requirements and optimally support you in your daily work.



PROGRESS NEEDS PIONEERS:

- / Product solutions for current and emerging trends
- / Innovative production technologies



DISCOVER
MORE:



We work hand in hand with you to
develop tailored solutions to suit
your requirements.



CUSTOMER FOCUS

We focus on you and
your goals

WE SUPPORT YOU SO WE CAN MAKE BIG THINGS HAPPEN TOGETHER

Everything we do and the way we do it is designed to make your work easier. By understanding your wishes and needs, as well as your requirements and workflows, we are able to offer you real added value and thus trust in long-term and sustainable customer relationships.

We are here to provide you with the best products and solutions

to suit your needs. Our dedicated team ensures that processes run smoothly and is always on hand when needed. Our multidisciplinary teams have many years of practical experience in science and laboratory work and understand the challenges of your daily work first hand. We are happy to use this expert knowledge to create training materials tailored to your specific needs.



WE OFFER:

- / Long-term customer relationship
- / Solutions to suit your requirements
- / Training and training materials



DISCOVER MORE:



We use over 60 years of experience
in the laboratory and healthcare
sectors, and our passion
for innovation, to meet your
requirements and grow with you.



EXPERTISE

In all relevant market trends

WE SPEAK YOUR LANGUAGE

We understand your tasks and their workflows with all their requirements. Our experts are at your side with more than 60 years of experience – from the development of products and solutions to support with regulatory issues or individual training.

More than 80% of our sales and application specialists have a background in natural sciences

and practical experience in the laboratory. We would be happy to advise you on site – at the place where your requirements arise. You can rely on our expertise to quickly integrate our products into your processes, because we know and understand your workflow. We are happy to support you in every step of your individual workflow, for example with product training, support documentation and service.



WHAT WE OFFER:

- / Comprehensive portfolio for your workflow needs
- / Experienced team of industry experts
- / Best practices for optimal efficiency



DISCOVER
MORE:



We help you to optimize your operations and processes so that you can always achieve the best possible results.



EFFICIENCY

We are your partner
for increasing
efficiency – inside
and outside the
laboratory

WE INCREASE YOUR EFFICIENCY

In the laboratory and day-to-day work, increasing efficiency is no longer just an option, but an essential requirement. As a strategic cooperation partner for day-to-day routine activities in and outside the laboratory, we support you in overcoming your current and future challenges.

We offer significant added value by supporting the efficien-

cy of your workflows. For example, our digital solutions create a connection between pre-bar-coded test tubes and patient data from the laboratory information system. This reliably rules out mix-ups and minimizes sampling errors. Your workflows can thus be measurably optimized, better results can be achieved with minimum time, resources, and work.



YOUR STRATEGIC PARTNER:

- / Increased efficiency
- / Digital solutions
- / Best practice recommendations



DISCOVER MORE:



We are committed to climate
protection and limiting global
warming to 1.5°C.



SUSTAINABILITY

The principle behind
our processes for
many years

BECAUSE IT'S REALLY IMPORTANT TO US

Sustainable action is the basic prerequisite for a future worth living in. As a plastics processing company, we strive to develop sustainable solutions that provide answers to the major challenges of our time.

Greiner AG's overarching "Blue-Plan" strategy focuses on the three pillars of climate, circular economy and people. Greiner


Bio-One is also working intensively on these sub-areas.

As a producer of plastic medical products, it is our declared long-term goal to combine ecological, economic and social aspects. We attach particular importance not only to your excellent supply of high-quality products, but also to energy-efficient and sustainable operation as a plastics producer.




OUR GOALS UNTIL 2030:

- / Comprehensive circular approach
- / Preparation of employees for the challenges of the future

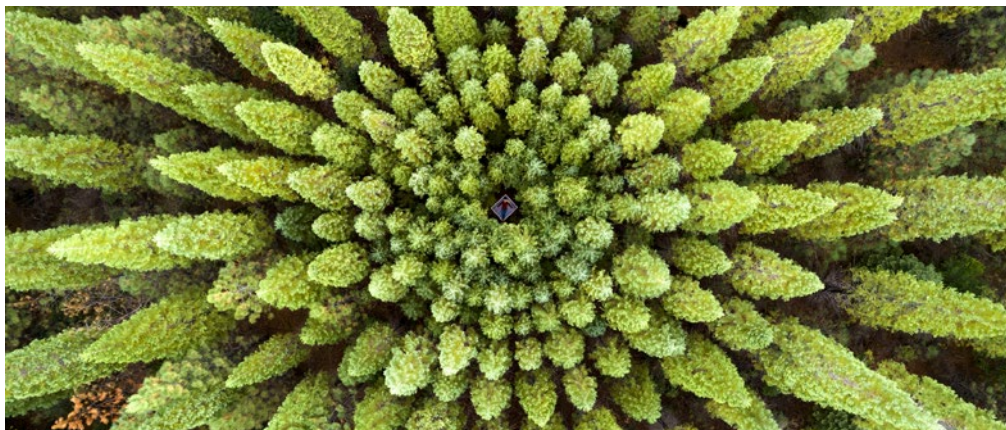

DISCOVER
MORE:

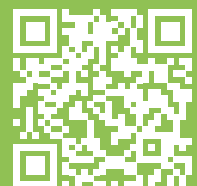



The innovative VACUETTE® Blood Collection Tubes made out of virtually unbreakable PET plastic have set the standard on today's market. They provide improved safety and hygiene during specimen collection procedures with the advantage of being as clear as glass.

VACUETTE®**BLOOD COLLECTION TUBES**

/ VACUETTE® Coagulation Tubes.....	40	/ VACUETTE® Crossmatch Tubes.....	63
9NC Coagulation Sodium Citrate 3.2 %	41	CAT Crossmatch Serum Clot Activator	64
9NC Coagulation Sodium Citrate 3.8 %	42	K3E Crossmatch K3EDTA	64
CTAD	42		
/ VACUETTE® Serum Tubes	43		
CAT Serum Clot Activator	44		
CAT Serum Separator Clot Activator	45		
CAT Serum Fast Separator	47		
/ VACUETTE® Heparin Tubes	48		
LH Lithium Heparin	49		
LH Lithium Heparin Separator	50		
NH Sodium Heparin	51		
/ VACUETTE® EDTA Tubes	52		
K2E K2EDTA	53		
K2E K2EDTA Separator	54		
K3E K3EDTA	54		
/ VACUETTE® No Additive Tubes.....	56		
Z No Additive	57		
/ VACUETTE® Glucose Tubes	58		
FX Sodium Fluoride / Potassium Oxalate	59		
FE Sodium Fluoride / K3EDTA	59		
FC Mix	60		
/ VACUETTE® Trace Element Tubes	61		
NH Trace Elements Sodium Heparin	62		

THE VACUUM SYSTEM: SIMPLE APPLICATION AND RELIABLE RESULTS

The main advantage of a user-friendly evacuated specimen collection system is that a clearly defined quantity of blood sample material can be collected without further handling requirements. VACUETTE® tubes can be changed quickly and hygienically during blood collection.

Most tubes contain chemical additives, which in combination with the pre-determined vacuum ensure the correct mixing ratio for the blood sample. The cap colors of the tubes are based on the internationally recommended color code so that different tube types can be recognized at first glance.





ADDITIONAL OPTICAL REC- OGNITION

thanks to different
colored rings on the
top of the cap.

All tubes with 13 mm diameter are available in two versions: as PREMIUM tubes with screw thread or as non-ridged tubes without screw thread. The VACUETTE® PREMIUM tube with screw thread allows hygienic opening and safe resealing of the tube during manual use. VACUETTE® non-ridged tubes can also be opened with a simple pulling motion.



PREMIUM Tubes with screw thread.

Due to the interaction of cap and thread on the PREMIUM tube, the VACUETTE® PREMIUM tube can be opened with a simple twisting movement with minimal effort.



Non-Ridged Tubes without screw thread.

Even without a thread, the non-ridged tubes are ideally opened with a twist and pull movement. All tubes with 16 mm diameter are offered with "non-ridged" pull-off caps as standard.



Of course, all tubes
have the CE-mark.

100% TRACK & TRACE WORLDWIDE

from production to
sample disposal

The item number is
encoded in a barcode
(GS1 DataMatrix)
together with
manufacturer-relevant
information such as
time and place of
production, lot number
or expiry date

Automated placement
of barcode labels
according to ISO
allows free view of
the sample material

Digital support in the
use of the correct sample
material during collection
possible through coded
sample type



Standardized
label and barcode
placement according
to ISO guideline
Code-128 for optimal
readability on all
common laboratory
equipment

High barcode quality
prevents smearing by
disinfectants or sweat

Support of GDPR
by pseudonymized
barcode = no patient
data on the tube



VACUETTE® barcode tubes are labeled during the manufacturing process.

PRE-BARCODED TUBES: THE NEW GENERATION FOR INCREASED SAFETY AND EFFICIENCY

Pre-barcode tubes enable a new digital sample collection process that improves the quality of laboratory results and enhances patient treatment. Ultimately, this helps increase efficiency and reduce costs.

The VACUETTE® barcode tubes form the basis for generating and processing information from the preanalytical process automatically using the appropriate software solution. The sample col-

lection process is made up of three simple steps: scanning the patient's wristband, scanning the tube and finally, handling in the laboratory.

Our pre-barcode tubes are in compliance with current data protection requirements and there is no need to apply extra labels, possibly blocking the viewing space, enabling optimal quality control.



**GUARANTEED
ISO BARCODE
QUALITY**

for automated
readability of
laboratory devices.



/ You can find more information about this product in our web-shop.



VACUETTE®

COAGULATION TUBES

VACUETTE® coagulation tubes with sodium citrate are used for determinations in coagulation tests.

Coagulation tubes from Greiner Bio-One are available in three versions: 9NC Coagulation Sodium Citrate 3.2 %, 9NC Coagulation Sodium Citrate 3.8 % and CTAD.

The additive trisodium citrate contained in the tube prevents

blood coagulation by binding Ca-ions.

Some of our VACUETTE® Coagulation Tubes are also known as "sandwich" tubes due to their double-wall-technology. The inner tube is made of PP (polypropylene) and prevents the evaporation of liquid. The outer tube is made of PET (polyethylene terephthalate) and ensures a stable vacuum for a correct draw volume.



STERILE

IVD

9NC Coagulation Sodium Citrate 3.2 %

VACUETTE® 9NC coagulation 3.2 % trisodium citrate (0.109 mol/l) tubes are filled with a buffered trisodium citrate solution. The choice of concentration depends on the laboratory's procedure. The mixing ratio is 1 part citrate to 9 parts blood. Selected versions are also available with a transparent label.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454321	2 ml	● light blue	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
474321	2 ml	● light blue	○ white	PREMIUM	13 x 75	transparent	no	50 / 1,200
484517	2 ml	● light blue	○ white	PREMIUM	13 x 75	paper	yes	50 / 1,200
454325	3 ml	● light blue	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
454327	3.5 ml	● light blue	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
474327	3.5 ml	● light blue	● black	PREMIUM	13 x 75	transparent	no	50 / 1,200
484518	3.5 ml	● light blue	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
486505	5 ml	● light blue	● black	PREMIUM	13 x 100	paper	yes	50 / 1,200
454322	2 ml	● light blue	○ white	non-ridged	13 x 75	paper	no	50 / 1,200
454334	3 ml	● light blue	● black	non-ridged	13 x 75	paper	no	50 / 1,200
484508	3 ml	● light blue	● black	non-ridged	13 x 75	paper	yes	50 / 1,200
454332	3.5 ml	● light blue	● black	non-ridged	13 x 75	paper	no	50 / 1,200



9NC Coagulation Sodium Citrate 3.8 %

VACUETTE® 9NC coagulation 3.8 % trisodium citrate (0.129 mol/l) tubes are filled with a buffered trisodium citrate solution. The choice of concentration depends on the laboratory's procedure. The mixing ratio is 1 part citrate to 9 parts blood.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454385	3 ml	● light blue	● black	PREMIUM	13 x 75	paper	no	50 / 1,200



CTAD

VACUETTE® CTAD tubes contain buffered citrate solutions, theophylline, adenosine and dipyridamole. As dipyridamole is extremely sensitive to light, the tubes are wrapped in white plastic film.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
474065	2 ml	● light blue	● yellow	PREMIUM	13 x 75	transparent	no	50 / 1,200
474064	3.5 ml	● light blue	● yellow	PREMIUM	13 x 75	transparent	no	50 / 1,200



/ You can find more information about this product in our web-shop.



VACUETTE® SERUM TUBES

VACUETTE® serum tubes are used for determinations in serum for routine clinical chemistry tests and hormones, serology, immunohematology and TDM.

Serum tubes from Greiner Bio-One are available in three versions: CAT Serum Clot Acti-

vator, CAT Serum Sep Clot Activator and CAT Serum Fast Separator.

The inner wall of each serum tube is coated with microscopic silica particles (clot activator main component SiO_2) which activate clotting when the tubes are gently inverted.



CAT Serum Clot Activator

VACUETTE® CAT serum clot activator tubes are coated with micronized silica particles that activate clotting when the tubes are gently inverted. Selected tubes are also available with a transparent label.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454098	1 ml	● red	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
454096	2 ml	● red	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
484515	2 ml	● red	○ white	PREMIUM	13 x 75	paper	yes	50 / 1,200
454092	4 ml	● red	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
474092	4 ml	● red	● black	PREMIUM	13 x 75	transparent	no	50 / 1,200
484525	4 ml	● red	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
456092	6 ml	● red	● black	PREMIUM	13 x 100	paper	no	50 / 1,200
476092	6 ml	● red	● black	PREMIUM	13 x 100	transparent	no	50 / 1,200
486504	6 ml	● red	● black	PREMIUM	13 x 100	paper	yes	50 / 1,200
454236	2 ml	● red	○ white	non-ridged	13 x 75	paper	no	50 / 1,200
454204	4 ml	● red	● black	non-ridged	13 x 75	paper	no	50 / 1,200
456089	6 ml	● red	● black	non-ridged	13 x 100	paper	no	50 / 1,200
456048	6 ml	● red	● black	non-ridged	13 x 100	transparent	no	50 / 1,200
455092	9 ml	● red	● black	non-ridged	16 x 100	paper	no	50 / 1,200



CAT Serum Separator Clot Activator

VACUETTE® CAT serum separator clot activator tubes contain a separating gel that migrates up to the interface between serum and clot during centrifugation, where it forms a stable barrier and separates the serum from the fibrin and cells.

Please note that depending on your country, different item numbers for the separator gel tubes may be available.

Item no.	Nominal volume	Cap color	Ring Color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454028	2.5 ml	red	yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454028A	2.5 ml	red	yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454067	3.5 ml	red	yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454067A	3.5 ml	red	yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454071	3.5 ml	red	yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454071A	3.5 ml	red	yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454078A	3.5 ml	red	yellow	non-ridged	13 x 75	transparent	no	50 / 1,200
454205A	2.5 ml	yellow	white	non-ridged	13 x 75	paper	no	50 / 1,200
454214	3.5 ml	gold	gold	PREMIUM	13 x 75	paper	no	50 / 1,200
454214A	3.5 ml	gold	gold	PREMIUM	13 x 75	paper	no	50 / 1,200
454225A	3.5 ml	red	yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454228	3.5 ml	gold	gold	non-ridged	13 x 75	paper	no	50 / 1,200
454228A	3.5 ml	gold	gold	non-ridged	13 x 75	paper	no	50 / 1,200
454243	2.5 ml	red	white	non-ridged	13 x 75	paper	no	50 / 1,200
454243A	2.5 ml	red	white	non-ridged	13 x 75	paper	no	50 / 1,200
454253A	3.5 ml	gold	white	non-ridged	13 x 75	paper	no	50 / 1,200
454260A	3.5 ml	white	yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454263A	3.5 ml	red	yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454319	3.5 ml	gold	black	non-ridged	13 x 75	paper	no	50 / 1,200
454319A	3.5 ml	gold	black	non-ridged	13 x 75	paper	no	50 / 1,200
455057	4 ml	red	yellow	non-ridged	16 x 100	paper	no	50 / 1,200
455057A	4 ml	red	yellow	non-ridged	16 x 100	paper	no	50 / 1,200
455071	8 ml	red	yellow	non-ridged	16 x 100	paper	no	50 / 1,200
455071A	8 ml	red	yellow	non-ridged	16 x 100	paper	no	50 / 1,200
456010	5 ml	gold	gold	PREMIUM	13 x 100	paper	no	50 / 1,200
456010A	5 ml	gold	gold	PREMIUM	13 x 100	paper	no	50 / 1,200
456069	5 ml	red	yellow	non-ridged	13 x 100	transparent	no	50 / 1,200
456069A	5 ml	red	yellow	non-ridged	13 x 100	transparent	no	50 / 1,200

Item no.	Nominal volume	Cap color	Ring Color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
456071	5 ml	● red	● yellow	PREMIUM	13 x 100	paper	no	50 / 1,200
456071A	5 ml	● red	● yellow	PREMIUM	13 x 100	paper	no	50 / 1,200
456073	5 ml	● red	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456073A	5 ml	● red	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456234	5 ml	● gold	● gold	non-ridged	13 x 100	transparent	no	50 / 1,200
456234A	5 ml	● gold	● gold	non-ridged	13 x 100	transparent	no	50 / 1,200
476010	5 ml	● gold	● gold	PREMIUM	13 x 100	transparent	no	50 / 1,200
476010A	5 ml	● gold	● gold	PREMIUM	13 x 100	transparent	no	50 / 1,200
476071	5 ml	● red	● yellow	PREMIUM	13 x 100	transparent	no	50 / 1,200
476071A	5 ml	● red	● yellow	PREMIUM	13 x 100	transparent	no	50 / 1,200
486503	5 ml	● red	● yellow	PREMIUM	13 x 100	paper	yes	50 / 1,200
474071A	3.5 ml	● red	● yellow	PREMIUM	13 x 75	transparent	no	50 / 1,200
474214	3.5 ml	● gold	● gold	PREMIUM	13 x 75	transparent	no	50 / 1,200
474214A	3.5 ml	● gold	● gold	PREMIUM	13 x 75	transparent	no	50 / 1,200
486500	5 ml	● gold	● gold	PREMIUM	13 x 100	paper	yes	50 / 1,200
454420A	3.5 ml	● gold	● gold	non-ridged	13 x 75	transparent	no	50 / 1,200
484513	2.5 ml	● red	● yellow	PREMIUM	13 x 75	paper	yes	50 / 1,200
484513A	2.5 ml	● red	● yellow	PREMIUM	13 x 75	paper	yes	50 / 1,200
484514	3.5 ml	● red	● yellow	PREMIUM	13 x 75	paper	yes	50 / 1,200
484514A	3.5 ml	● red	● yellow	PREMIUM	13 x 75	paper	yes	50 / 1,200
485502	8 ml	● red	● yellow	non-ridged	16 x 100	paper	yes	50 / 1,200
455072	8 ml	● gold	● gold	non-ridged	16 x 100	paper	no	50 / 1,200
455072A	8 ml	● gold	● gold	non-ridged	16 x 100	paper	no	50 / 1,200
455099A	7 ml	● red	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200
456004A	5 ml	○ white	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456005A	5 ml	● yellow	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456006A	5 ml	● royal blue	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456018	5 ml	● gold	● gold	non-ridged	13 x 100	paper	no	50 / 1,200
456018A	5 ml	● gold	● gold	non-ridged	13 x 100	paper	no	50 / 1,200
456019A	5 ml	● pink	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456214A	3.5 ml	● red	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
485514	4 ml	● red	● yellow	non-ridged	16 x 100	paper	yes	50 / 1,200
485514A	4 ml	● red	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200
486004A	5 ml	○ white	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
486005A	5 ml	● red	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
486006A	5 ml	● red	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
486007A	5 ml	● yellow	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200



STERILE

IVD

CAT Serum Fast Separator

VACUETTE® CAT serum fast separator tubes contain a separating gel at the bottom of the tube. In contrast to conventional serum tubes, the Serum Fast Tube contains additional thrombin to accelerate the clotting process.

Please note that depending on your country, different item numbers for the separator gel tubes may be available.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454592	3.5 ml	orange	yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
456309	5 ml	orange	yellow	PREMIUM	13 x 100	paper	no	50 / 1,200
486509	5 ml	orange	yellow	PREMIUM	13 x 100	paper	yes	50 / 1,200
454593	3.5 ml	orange	yellow	non-ridged	13 x 75	paper	no	50 / 1,200
456313	5 ml	orange	yellow	non-ridged	13 x 100	paper	no	50 / 1,200
454591	3.5 ml	orange	yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454591A	3.5 ml	orange	yellow	non-ridged	13 x 75	paper	no	50 / 1,200
456310A	5 ml	orange	yellow	non-ridged	13 x 100	paper	no	50 / 1,200



/ You can find more information about this product in our web-shop.



VACUETTE® HEPARIN TUBES

VACUETTE® heparin tubes are used for plasma determinations of routine clinical chemistry tests.

Heparin tubes from Greiner Bio-One are available in three different versions: LH Lithium

Heparin, LH Lithium Heparin Sep and NH Sodium Heparin.

The heparin additive activates ATIII, blocking the coagulation cascade, resulting in whole blood / plasma sample.



STERILE

IVD

LH Lithium Heparin

VACUETTE® LH lithium heparin tubes are coated with lithium heparin on the inside. The anticoagulant heparin activates antithrombin, which blocks the clotting cascade and thus produces a whole blood/plasma sample. Selected tubes are also available with a transparent label.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454081	1 ml	● green	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
454089	2 ml	● green	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
484589	2 ml	● green	○ white	PREMIUM	13 x 75	paper	yes	50 / 1,200
474084	4 ml	● green	● black	PREMIUM	13 x 75	transparent	no	50 / 1,200
484522	4 ml	● green	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
486508	6 ml	● green	● black	PREMIUM	13 x 100	paper	yes	50 / 1,200
456084	6 ml	● green	● black	PREMIUM	13 x 100	paper	no	50 / 1,200
454237	2 ml	● green	○ white	non-ridged	13 x 75	paper	no	50 / 1,200
455222	4 ml	● green	● green	non-ridged	16 x 100	paper	no	50 / 1,200
485222	4 ml	● green	● black	non-ridged	16 x 100	paper	yes	50 / 1,200
456088	6 ml	● green	● black	non-ridged	13 x 100	paper	no	50 / 1,200
455084	9 ml	● green	● black	non-ridged	16 x 100	paper	no	50 / 1,200



LH Lithium Heparin Separator

VACUETTE® LH lithium heparin separator tubes also contain a separating gel that moves upwards during centrifugation, where it then creates a stable barrier that separates the plasma from the cells. Selected tubes are also available with a transparent label.

Please note that depending on your country, different item numbers for the separator gel tubes may be available.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454008	3.5 ml	● green	● yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454008A	3.5 ml	● green	● yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454046A	2.5 ml	● green	● yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454083	3.5 ml	● green	● yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454083A	3.5 ml	● green	● yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
454483A	3.5 ml	● green	● yellow	non-ridged	13 x 75	paper	no	50 / 1,200
454489	3.5 ml	● green	● yellow	non-ridged	13 x 75	transparent	no	50 / 1,200
454489A	3.5 ml	● mint green	● yellow	non-ridged	13 x 75	transparent	no	50 / 1,200
454491A	3 ml	● green	● yellow	non-ridged	13 x 75	transparent	no	50 / 1,200
454493A	3 ml	● mint green	● yellow	non-ridged	13 x 75	transparent	no	50 / 1,200
455083	8 ml	● green	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200
455083A	8 ml	● green	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200
455229	4 ml	● green	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200
455229A	4 ml	● green	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200
456083	5 ml	● green	● yellow	PREMIUM	13 x 100	paper	no	50 / 1,200
456083A	5 ml	● green	● yellow	PREMIUM	13 x 100	paper	no	50 / 1,200
456087	5 ml	● green	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456087A	5 ml	● green	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456216	5 ml	● green	● yellow	non-ridged	13 x 100	transparent	no	50 / 1,200
456216A	5 ml	● green	● yellow	non-ridged	13 x 100	transparent	no	50 / 1,200
456287A	3.5 ml	● green	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
456305A	5 ml	● mint green	● yellow	non-ridged	13 x 100	transparent	no	50 / 1,200
474083A	3.5 ml	● green	● yellow	PREMIUM	13 x 75	transparent	no	50 / 1,200
484523	3.5 ml	● green	● yellow	non-ridged	13 x 75	paper	yes	50 / 1,200
484523A	3.5 ml	● green	● yellow	PREMIUM	13 x 75	paper	yes	50 / 1,200
476083	5 ml	● green	● yellow	PREMIUM	13 x 100	transparent	no	50 / 1,200
485513	4 ml	● green	● yellow	non-ridged	16 x 100	paper	yes	50 / 1,200
486507	5 ml	● green	● yellow	non-ridged	13 x 100	paper	yes	50 / 1,200

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454247A	3 ml	● green	● yellow	non-ridged	13 x 75	paper	no	50 / 1,200
484030A	3.5 ml	● green	● red	non-ridged	13 x 75	paper	yes	50 / 1,200
484031A	3.5 ml	● green	● yellow	non-ridged	13 x 75	paper	yes	50 / 1,200



NH Sodium Heparin

The inner wall of the VACUETTE® NH sodium heparin tubes is coated with spray-dried sodium heparin. The anticoagulant heparin activates antithrombin, which blocks the clotting cascade, producing a whole blood/ plasma sample.



Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
484531	4 ml	● green	● green	PREMIUM	13 x 75	paper	yes	50 / 1,200
455051	9 ml	● green	● green	non-ridged	16 x 100	paper	no	50 / 1,200



/ You can find more information about this product in our web-shop.



VACUETTE® EDTA TUBES

VACUETTE® EDTA tubes are used for determinations in EDTA whole blood for hematology and routine immunohematology. K2E K2EDTA Sep tubes are also available for molecular diagnostics and viral load detection.

EDTA tubes from Greiner Bio-One are available in three versions: K2E K2EDTA, K2E K2EDTA Sep-

arator and K3E K3 EDTA. The additive EDTA (ethylenediaminetetraacetic acid) contained in the tube interferes with calcium and thus blockades the coagulation cascade.

- / Di-Potassium K2EDTA
- / Tri-Potassium K3EDTA



K2E K2EDTA

VACUETTE® K2E K2EDTA tubes contain di-potassium K2EDTA which blocks the coagulation cascade. Selected tubes are also available with transparent labels.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
484586	2 ml	● lavender	○ white	PREMIUM	13 x 75	paper	yes	50 / 1,200
454024	2 ml	● lavender	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
454020	3 ml	● lavender	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
484526	3 ml	● lavender	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
454023	4 ml	● lavender	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
456023	6 ml	● lavender	● black	PREMIUM	13 x 100	paper	no	50 / 1,200
454047	2 ml	● lavender	○ white	non-ridged	13 x 75	paper	no	50 / 1,200
454246	3 ml	● lavender	● black	non-ridged	13 x 75	paper	no	50 / 1,200
454411	3 ml	● lavender	● black	non-ridged	13 x 75	transparent	no	50 / 1,200
454209	4 ml	● lavender	● black	non-ridged	13 x 75	paper	no	50 / 1,200
454410	4 ml	● lavender	● black	non-ridged	13 x 75	transparent	no	50 / 1,200
484507	4 ml	● lavender	● black	non-ridged	13 x 75	paper	yes	50 / 1,200
456043	6 ml	● lavender	● black	non-ridged	13 x 100	paper	no	50 / 1,200
455045	9 ml	● lavender	● black	non-ridged	16 x 100	paper	no	50 / 1,200



K2E K2EDTA Separator

VACUETTE® K2E K2EDTA separator tubes contain an additional separating gel and are used for plasma testing in molecular diagnostics and for virus load determination. Selected tubes are also available with a transparent label.

Item no.	Nominal volume	Cap color	Ring Color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454235	3.5 ml	● lavender	● yellow	PREMIUM	13 x 75	paper	no	50 / 1,200
456211	5 ml	○ white	● yellow	PREMIUM	13 x 100	paper	no	50 / 1,200
456058	5 ml	○ white	● yellow	non-ridged	13 x 100	paper	no	50 / 1,200
455040	8 ml	● lavender	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200
455021	8 ml	○ white	● yellow	non-ridged	16 x 100	paper	no	50 / 1,200



K3E K3EDTA

VACUETTE® K3E K3EDTA tubes contain tri-potassium K3EDTA which blocks the coagulation cascade. Selected tubes are also available with transparent labels.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
484509	1 ml	● lavender	○ white	PREMIUM	13 x 75	paper	yes	50 / 1,200
454087	2 ml	● lavender	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
474087	2 ml	● lavender	○ white	PREMIUM	13 x 75	transparent	no	50 / 1,200
484512	2 ml	● lavender	○ white	PREMIUM	13 x 75	paper	yes	50 / 1,200

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454086	3 ml	● lavender	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
474086	3 ml	● lavender	● black	PREMIUM	13 x 75	transparent	no	50 / 1,200
484511	3 ml	● lavender	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
454036	4 ml	● lavender	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
474036	4 ml	● lavender	● black	PREMIUM	13 x 75	transparent	no	50 / 1,200
484510	4 ml	● lavender	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
456036	6 ml	● lavender	● black	PREMIUM	13 x 100	paper	no	50 / 1,200
486502	6 ml	● lavender	● black	PREMIUM	13 x 100	paper	yes	50 / 1,200
454222	2 ml	● lavender	○ white	non-ridged	13 x 75	paper	no	50 / 1,200
454217	3 ml	● lavender	● black	non-ridged	13 x 75	paper	no	50 / 1,200
454287	3 ml	● lavender	● black	non-ridged	13 x 75	transparent	no	50 / 1,200
454021	4 ml	● lavender	● black	non-ridged	13 x 75	paper	no	50 / 1,200
456038	6 ml	● lavender	● black	non-ridged	13 x 100	paper	no	50 / 1,200
455036	9 ml	● lavender	● black	non-ridged	16 x 100	paper	no	50 / 1,200
485501	9 ml	● lavender	● black	non-ridged	16 x 100	paper	yes	50 / 1,200



/ You can find more information about this product in our web-shop.



VACUETTE®

NO ADDITIVE TUBES

VACUETTE® no additive tubes are used as discard tubes or for the collection of blood.

The tubes without additive from Greiner Bio-One are available in

the version Z no additive. There is no chemical additive in the tube.

**Z No Additive**

VACUETTE® tubes without additives contain no chemical additives, but are evacuated and sterile on the inside.



Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454088	2 ml	○ white	○ white	PREMIUM	13 x 75	paper	no	50 / 1,200
454001	4 ml	○ white	● black	PREMIUM	13 x 75	paper	no	50 / 1,200
484530	4 ml	○ white	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
456001	6 ml	○ white	● black	PREMIUM	13 x 100	paper	no	50 / 1,200
454241	3 ml	○ white	● black	non-ridged	13 x 75	paper	no	50 / 1,200
456085	6 ml	○ white	● black	non-ridged	13 x 100	paper	no	50 / 1,200
455001	9 ml	○ white	● black	non-ridged	16 x 100	paper	no	50 / 1,200
485506	9 ml	○ white	● black	non-ridged	16 x 100	paper	yes	50 / 1,200



/ You can find more information about this product in our web-shop.



VACUETTE®

GLYCOLYTIC INHIBITOR TUBES

VACUETTE® sodium fluoride/K3EDTA or sodium fluoride/potassium oxalate tubes are suitable for the analysis of glucose concentration, including the determination of lactate.

The VACUETTE® FX Mix tubes are very effective for immediate stabilization of the in vivo glucose concentration.

The powder additive in the FC mix tube lowers the pH and blocks the pH-dependent enzymes that would be active in the initial phase of the glycolysis cascade. The combination of the citrate-fluoride-EDTA additive in the tube stabilizes the in vivo concentration of glucose in the sample from the beginning.



STERILE

IVD

FX Sodium Fluoride / Potassium Oxalate

These VACUETTE® tubes contain sodium fluoride as a stabilizer and the additive potassium oxalate as an anticoagulant.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
484528	2 ml	● grey	○ white	non-ridged	13 x 75	paper	yes	50 / 1,200
454297	4 ml	● grey	● black	non-ridged	13 x 75	paper	no	50 / 1,200



STERILE

IVD

FE Sodium Fluoride / K3EDTA

In addition to sodium fluoride as a stabilizer, the VACUETTE® tubes also contain the additive K3EDTA as an anticoagulant.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
484520	2 ml	● grey	○ white	PREMIUM	13 x 75	paper	yes	50 / 1,200
484527	4 ml	● grey	● black	PREMIUM	13 x 75	paper	yes	50 / 1,200
454221	2 ml	● grey	○ white	non-ridged	13 x 75	paper	no	50 / 1,200

VACUETTE® Blood Collection Tubes

VACUETTE® Glycolytic Inhibitor Tubes



STERILE

IVD

FC Mix

FC Mix tubes contain buffered Na₂ EDTA, sodium fluoride, citric acid and sodium citrate. This additive mixture reduces the pH value and blocks the pH-dependent enzymes that would be active in the initial stage of the glycolysis cascade.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
454510	2 ml	● pink	○ white	non-ridged	13 x 75	transparent	no	50 / 1,200
454511	2 ml	● grey	○ white	non-ridged	13 x 75	transparent	no	50 / 1,200
454513	3 ml	● pink	● black	non-ridged	13 x 75	transparent	no	50 / 1,200
454514	3 ml	● grey	● black	non-ridged	13 x 75	transparent	no	50 / 1,200



/ You can find more information about this product in our web-shop.



VACUETTE® TRACE ELEMENT TUBES

VACUETTE® trace element tubes are used to determine trace elements such as arsenic, copper, mercury and many more.

The trace element tube from Greiner Bio-One is available in the version NH trace elements sodium heparin.

Trace element tubes contain sodium heparin and are used to investigate trace elements. Prior to trace element determination, all equipment used during collection, transport and storage should be checked. A blank measurement must first be carried out for each tube batch.



NH Trace Elements Sodium Heparin

The tubes are suitable for testing most trace elements. NH trace elements sodium heparin tubes contain the anti-coagulant sodium heparin.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
456080	6 ml	● royal blue	● black	PREMIUM	13 x 100	paper	no	50 / 1,200
456275	6 ml	● royal blue	● black	non-ridged	13 x 100	paper	no	50 / 1,200



/ You can find more information about this product in our web-shop.



VACUETTE® CROSSMATCH TUBES

VACUETTE® crossmatch tubes are used for serological compatibility testing.

Two versions are available. CAT serum crossmatch contains clot activator for crossmatching with serum samples and the and K3E

EDTA is used for crossmatching with whole blood.



CAT Crossmatch Serum Clot Activator

CAT serum clot activator crossmatch tubes contain a clot activator to enable crossmatching in serum.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
456097	6 ml	● pink	● black	non-ridged	13 x 100	paper	no	50 / 1,200



K3E Crossmatch K3EDTA

K3E K3EDTA crossmatch tubes contain tri-potassium K3EDTA to enable crossmatching in whole blood.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Label	Barcode	Qty. inner / outer
456093	6 ml	● pink	● black	PREMIUM	13 x 100	paper	no	50 / 1,200
456052	6 ml	● pink	● black	non-ridged	13 x 100	paper	no	50 / 1,200




MiniCollect® offers a gentle way to collect small blood samples for a wide range of analyses. The system is perfectly suited for young children, geriatric patients as well as patients with fragile veins. It is also recommended in situations where venous blood collection proves particularly difficult.

MINICOLLECT®**BLOOD COLLECTION TUBES**

/ MiniCollect® Coagulation Tubes	72
9NC Coagulation Sodium Citrate 3.2 %	73
/ MiniCollect® Serum Tubes	74
CAT Serum Clot Activator	75
CAT Serum Separator Clot Activator	75
/ MiniCollect® Heparin Tubes	76
LH Lithium Heparin	77
LH Lithium Heparin Separator	77
/ MiniCollect® EDTA Tubes	78
K2E K2EDTA	79
K3E K3EDTA	79
/ MiniCollect® Glycolytic Inhibitor Tubes..	80
FX Sodium Fluoride / Potassium Oxalate	81

for collecting, transporting & processing capillary blood

for collecting, transporting & processing capillary blood

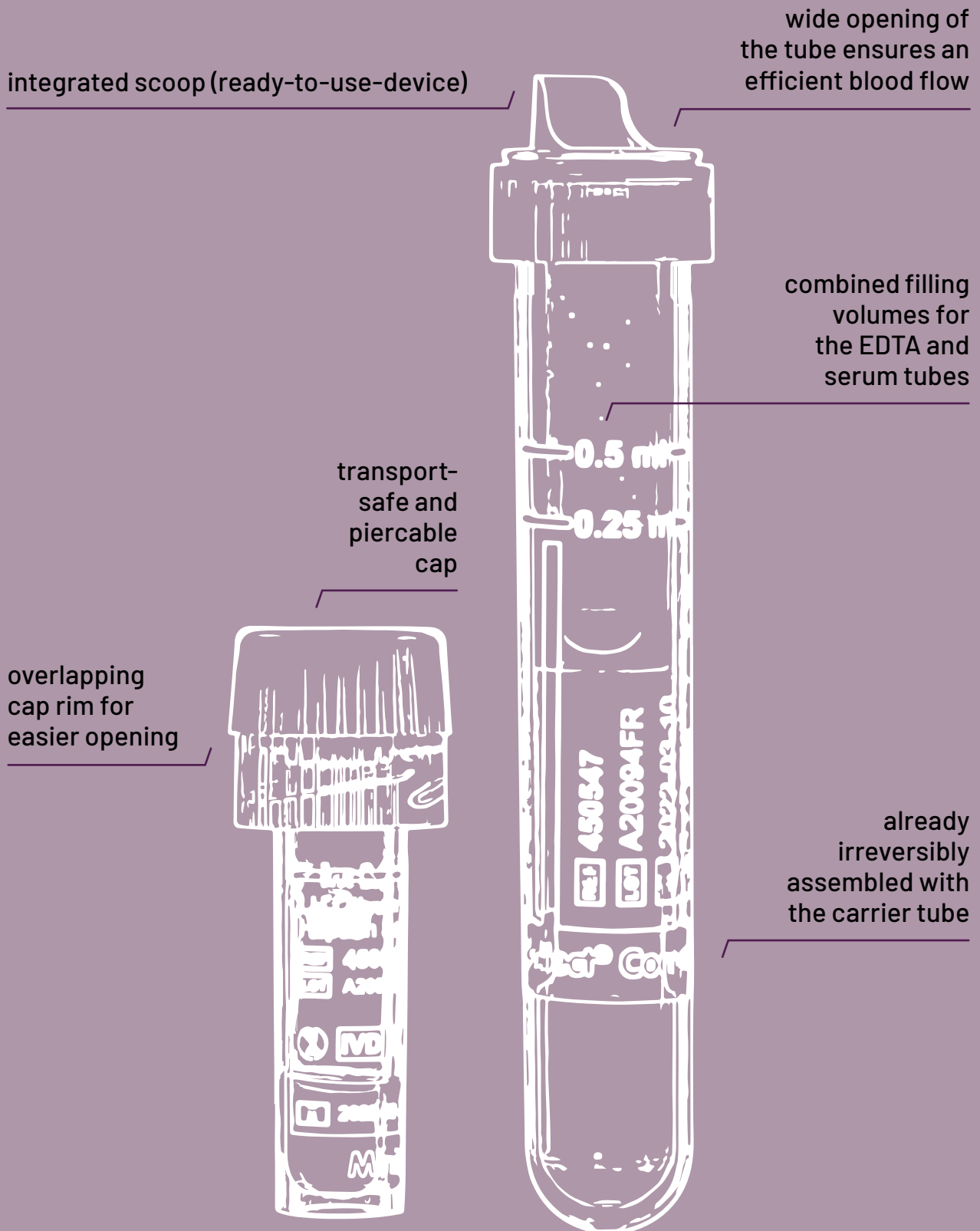



Easy and hygienic sample collection with MiniCollect®.

A GENTLE WAY TO COLLECT SMALL BLOOD SAMPLES

MiniCollect® offers a gentle way to collect small blood samples for a wide range of analyses. The system is perfectly suited for young children, geriatric patients as well as patients with fragile veins.

MiniCollect® is also recommended in situations where venous blood collection proves particularly difficult or contraindicated. For delicate vein con-

ditions, patients who have to take blood samples very often or people with severe burns, capillary blood collection with the MiniCollect® system offers a great alternative to venous blood collection.

The MiniCollect® tube includes an integrated scoop, allowing droplets of blood to be collected easily and hygienically. The additional use of a capillary or funnel is not necessary anymore.

MAXIMUM TUBE OPENING FOR MINIMUM ADHESION.

- / wide opening of the tube ensures an efficient blood flow
- / increased diameter facilitates the mixing of the sample
- / minimizes adhesion to the wall of the tube



MINICOLLECT® WITH CARRIER TUBE

for your automated laboratory process

MINIMAL BLOOD VOLUME YET MAXIMUM SPECTRUM OF ANALYSIS

The capillary blood sample is becoming increasingly popular as specimen material due to its simpler and less invasive sampling procedure. One of the advantages is the low blood volumes.

Small quantities are sufficient to enable a variety of parameters to be tested. Capillary blood collec-

tion supports the Patient Blood Management approach: As little blood as possible, as much as necessary. If for some reason venous blood collection is not possible, capillary blood collection could be a promising alternative. For children and patients who require frequent sampling it is essential that the volume of blood collected is minimized.





TWO VOLUMES ONE TUBE

Combined filling
volumes simplify the
sample processing.

Combined filling volumes for the EDTA and serum tubes simplify sample processing and reduce logistic demands. To ensure a correct mixing ratio for the capillary blood sample with the additive, the tubes just need to be filled to between the fill marks.



Transparent labels and color-coded caps.

The various colored caps reflect the international color codes according to ISO 6710. The colored label text ensures that the type of tube can be identified even when the cap is removed. The transparent label allows a visual inspection of the sample material.



Transport-safe caps.

The caps are virtually leakproof and can be transported via pneumatic post without loss of sample material. MiniCollect® tubes are leakproof according to ISO 6710 Annex C.



Ready to use device.

Due to the integrated scoop, no additional capillary or funnel is needed for sample collection. This increases efficiency compared to conventional capillary blood collection systems. Inventory management and planning can be completed easily and efficiently as no additional accessories are required.



/ You can find more information about this product in our web-shop.



MINICOLLECT® COAGULATION TUBES

MiniCollect® Coagulation Tubes with sodium citrate are used for determinations in coagulation tests.

Coagulation Tubes from Greiner

Bio-One are available in one version with sodium citrate 3.2%. The additive trisodium citrate contained in the tube prevents blood coagulation by binding Ca-ions.

**9NC Coagulation Sodium Citrate 3.2 %**

MiniCollect® 9NC Coagulation 3.2 % Trisodium Citrate Tubes are filled with a buffered trisodium citrate solution at a concentration of 0.109 mol/l and are used to collect venous blood samples anticoagulated with citrate for coagulation tests.

Model: MiniCollect® Tube, Nominal volume: 1 ml

Item no.	Cap color	Label	Barcode	Qty. inner / outer
450539	● light blue	transparent	no	500 / 2,000
480539	● light blue	paper	yes	500 / 2,000



/ You can find more information about this product in our web-shop.



MINICOLLECT® SERUM TUBES

MiniCollect® Serum Tubes are used for determinations in clinical chemistry from serum.

Serum tubes from Greiner Bio-One are available in two versions: CAT Serum (without separating gel) and CAT Serum

Separator (with separating gel). The inner wall of each serum tube is coated with microscopic silica particles (clot activator main component SiO_2) which activate clotting when the tubes are gently inverted.



CAT Serum Clot Activator

MiniCollect® CAT Serum Clot Activator tubes are available with 0.5 to 1 ml fill volume. The serum tubes are available as standard MiniCollect® tubes or as MiniCollect® Complete tubes in pre-assembled carrier tubes.

Volume range: 0.5 - 1 ml

Item no.	Model	Cap color	Label	Barcode	Qty. inner / outer
450534	MiniCollect® Tube	red	transparent	no	500 / 2,000
450549	MiniCollect® Complete	red	transparent	no	500 / 2,000
480548	MiniCollect® Complete	red	paper	yes	500 / 2,000



CAT Serum Separator Clot Activator

MiniCollect® CAT Serum Separator Clot Activator Tubes are available with 0.5 to 0.8 ml fill volume and contain an additional separating gel. The serum tubes are available as standard MiniCollect® tubes or as MiniCollect® Complete tubes in pre-assembled carrier tubes.

Volume range: 0.5 - 0.8 ml

Item no.	Model	Cap color	Label	Barcode	Qty. inner / outer
450533	MiniCollect® Tube	gold	transparent	no	500 / 2,000
450548	MiniCollect® Complete	gold	transparent	no	500 / 2,000
450568	MiniCollect® Complete	gold	paper	no	500 / 2,000



/ You can find more information about this product in our web-shop.



MINICOLLECT® HEPARIN TUBES

MiniCollect® Heparin Tubes are used for determinations in clinical chemistry from plasma.

Greiner Bio-One lithium heparin tubes are available in two versions: Lithium Heparin (without separating gel), Lithium Heparin

Separator (with separating gel).

The heparin additive activates ATIII, blocking the coagulation cascade, resulting in whole blood / plasma sample.



LH Lithium Heparin

MiniCollect® LH Lithium Heparin tubes are coated with lithium heparin on the inside. The anticoagulant heparin activates antithrombin, thus blocking the coagulation cascade and producing a whole blood/plasma sample. This is ideal for rapid analyses and blood analyses of patients on anticoagulant therapy.

Item no.	Model	Nominal volume	Cap color	Label	Barcode	Qty. inner / outer
450536	MiniCollect® Tube	0.5 ml	● green	transparent	no	500 / 2,000
450537	MiniCollect® Tube	1 ml	● green	transparent	no	500 / 2,000
480550	MiniCollect® Complete	0.5 ml	● green	paper	yes	500 / 2,000
450551	MiniCollect® Complete	1 ml	● green	transparent	no	500 / 2,000



LH Lithium Heparin Separator

MiniCollect® LH Lithium Heparin Separator Tubes also contain a separating gel. During centrifugation, the separating gel migrates upwards and acts as a barrier between the cells and the plasma.

Nominal volume: 0.8 ml

Item no.	Model	Cap color	Label	Barcode	Qty. inner / outer
450535	MiniCollect® Tube	● light green	transparent	no	500 / 2,000
450550	MiniCollect® Complete	● light green	transparent	no	500 / 2,000
480560	MiniCollect® Complete	● light green	paper	yes	500 / 2,000



/ You can find more information about this product in our web-shop.



MINICOLLECT® EDTA TUBES

MiniCollect® EDTA Tubes are used for the determination of the following hematological parameters: WBC, RBC, HGB, HCT, MCV, MCH, MCHC, platelets, RDW, lymphocytes, neutrophils, monocytes, eosinophils and basophils.

EDTA Tubes from Greiner Bio-One are available in two versions: K3EDTA and K2EDTA. The additive dipotassium EDTA (ethylenediaminetetraacetic acid) contained in the tube binds with calcium ions and thus blocks the coagulation cascade.

**K2E K2EDTA**

MiniCollect® K2E K2EDTA tubes contain di-potassium K2EDTA and are available with 0.25 to 0.5 ml fill volume. K2EDTA tubes are available as standard MiniCollect® tubes or as MiniCollect® Complete tubes in pre-assembled carrier tubes.

Volume range: 0.25 - 0.5 ml

Item no.	Model	Cap color	Label	Barcode	Qty. inner / outer
450532	MiniCollect® Tube	● lavender	transparent	no	500 / 2,000
450547	MiniCollect® Complete	● lavender	transparent	no	500 / 2,000
480547	MiniCollect® Complete	● lavender	paper	yes	500 / 2,000

**K3E K3EDTA**

MiniCollect® K3E K3EDTA tubes contain tri-potassium K3EDTA and are available in 0.25 to 0.5 ml and 1 ml volumes. K3EDTA tubes are available as standard MiniCollect® tubes or as MiniCollect® Complete tubes in pre-assembled carrier tubes.

Item no.	Model	Nominal volume	Volume range	Cap color	Label	Barcode	Qty. inner / outer
450530	MiniCollect® Tube		0.25 - 0.5 ml	● lavender	transparent	no	500 / 2,000
450531	MiniCollect® Tube	1 ml		● lavender	transparent	no	500 / 2,000
450545	MiniCollect® Complete		0.25 - 0.5 ml	● lavender	transparent	no	500 / 2,000
480546	MiniCollect® Complete		0.25 - 0.5 ml	● lavender	paper	yes	500 / 2,000
450546	MiniCollect® Complete	1 ml		● lavender	transparent	no	500 / 2,000



/ You can find more information about this product in our web-shop.



MINICOLLECT®

GLYCOLYTIC INHIBITOR TUBES

MiniCollect® Glycolytic Inhibitor Tubes are used to determine glucose and lactate.

Glycolytic inhibitor tubes from Greiner Bio-One contain the ad-

ditive sodium fluoride/potassium oxalate.

The spray-dried additive in the MiniCollect® Glycolytic Inhibitor Tubes ensures easy mixing of blood and additive.



FX Sodium Fluoride / Potassium Oxalate

These MiniCollect® tubes are available with 0.25 ml and 0.5 ml filling volume. The glycolytic inhibitor tubes are available as standard MiniCollect® tubes or as MiniCollect® Complete tubes in pre-assembled carrier tubes. The Complete version is only available with 0.25 ml filling volume.

Item no.	Model	Nominal volume	Cap color	Label	Barcode	Qty. inner / outer
450540	MiniCollect® Tube	0.25 ml	● grey	transparent	no	500 / 2,000
450541	MiniCollect® Tube	0.5 ml	● grey	transparent	no	500 / 2,000
450552	MiniCollect® Complete	0.25 ml	● grey	transparent	no	500 / 2,000
480549	MiniCollect® Complete	0.25 ml	● grey	paper	yes	500 / 2,000



One of the most effective ways to protect yourself against needlestick injuries is to use safety products in your daily work. We offer a wide range of safety products so that the most suitable products can be selected for each market and customer.


SAFETY PRODUCTS

/ VACUETTE® QUICKSHIELD	88
VACUETTE® QUICKSHIELD Safety Tube Holder	89
VACUETTE® QUICKSHIELD Complete	89
VACUETTE® QUICKSHIELD Complete PLUS	90
 / VACUETTE® EVOPROTECT	91
VACUETTE® EVOPROTECT SAFETY Blood Collection/Infusion Set	92
VACUETTE® EVOPROTECT SAFETY Blood Collection Set + Luer Adapter	92
VACUETTE® EVOPROTECT SAFETY Blood Collection Set + Holder	93
 / SAFETY Blood Collection Set	94
SAFETY Blood Collection/Infusion Set	95
SAFETY Blood Collection Set + Luer Adapter ..	95
SAFETY Blood Collection Set + Holder	96
SAFETY Blood Collection Set + Blood Culture Holder	96
 / MiniCollect® Safety Lancets	97
MiniCollect® Lancelino Safety Lancets with Blade	98
MiniCollect® Lancelino Safety Lancets with Needle	98
MiniCollect® Safety Lancets with Blade	99
MiniCollect® Safety Lancets with Needle	99

AN END TO NEEDLESTICK INJURIES

A large number of all occupationally acquired infectious diseases in the health sector are due to needlestick injuries. Needlestick injuries pose by far the highest risk of infection. More than 50 % of all registered cases affect nursing staff, followed by doctors and laboratory staff.

However, this does not mean that needlestick injuries represent an unavoidable occupational risk. The correct use of VACUETTE® safety products can help reduce this risk.

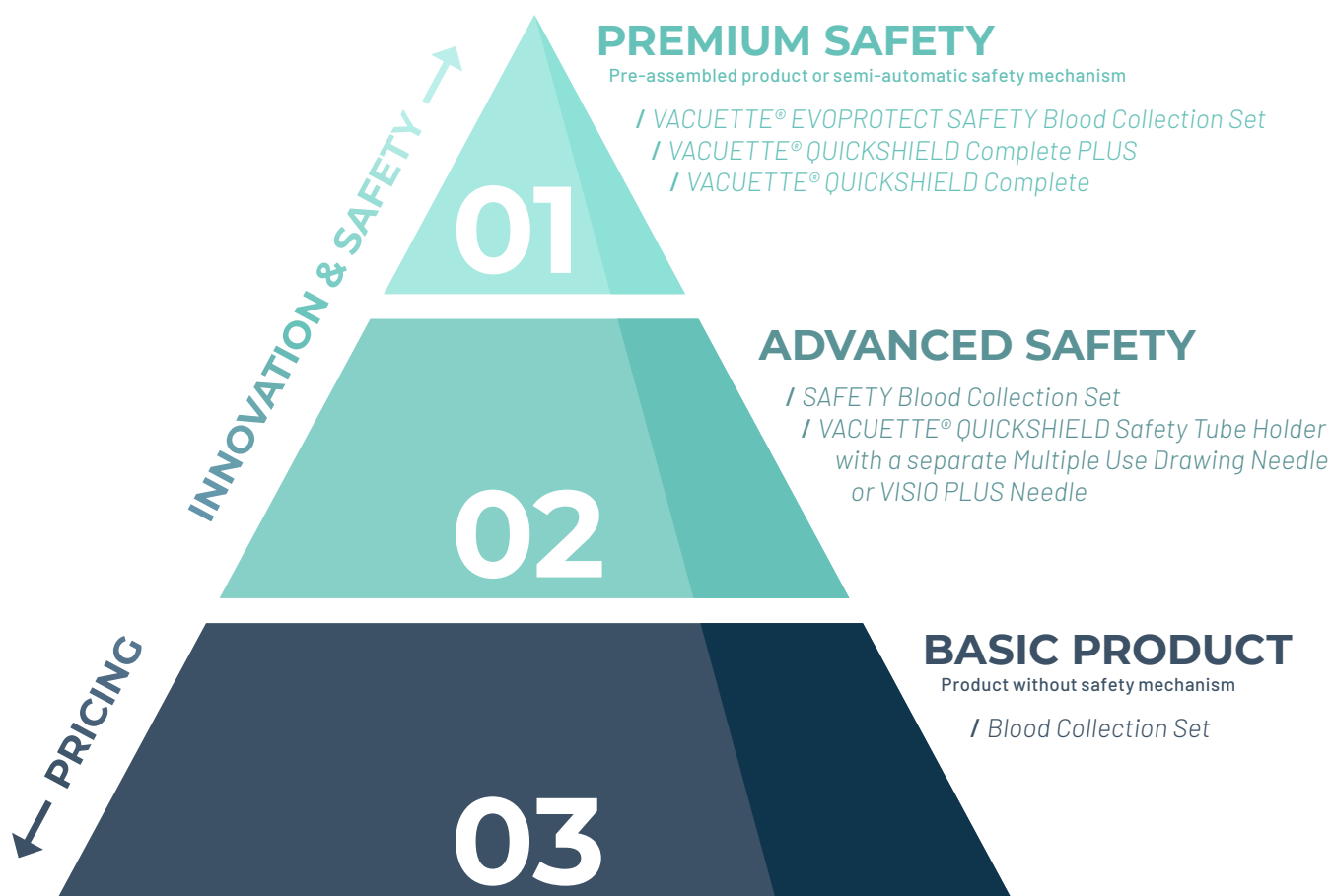



In our range you will find the right product for every safety level and your individual needs.

From products without a safety mechanism to premium safety products, we offer you the product for your individual needs. On the following pages you will get an overview of our products to protect yourself in the best possible way against needlestick injuries.



Standard venipuncture.

For daily blood collection we recommend the use of our VACUETTE® QUICKSHIELD Complete PLUS.

*VACUETTE® QUICKSHIELD
Complete PLUS*



Difficult vein conditions.

All our safety blood collection sets are also suitable for patients with difficult venous conditions.

*VACUETTE® EVOPROTECT or
SAFETY Blood Collection Set*



Blood cultures.

For the collection of blood cultures we can recommend our SAFETY Blood Collection Set with pre-assembled blood culture holder.

*SAFETY Blood Collection Set
with Blood Culture Holder
(21G or 23G)*

THREE TYPES OF LANCETS

for both young and
elderly patients

available with
needle or blade

ergonomic
design
for easy
handling


button
activated
safety lancet



contact
activated
safety lancet

available with
needle or blade

ergonomic
design for
easy handling

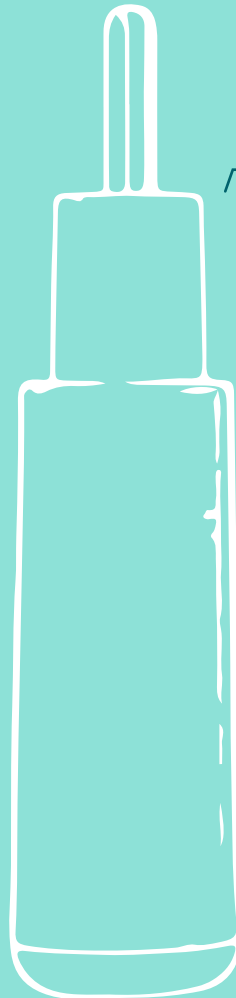



The blade/needle is automatically retracted after the puncture and is safely enclosed within the plastic housing.

WIDE RANGE OF SAFETY LANCETS

The MiniCollect® lancets are available with a wide range of puncture depths and blade sizes / needle gauges to ensure that the puncture wound can be kept as small as possible.

The safety mechanism ensures that the blade/needle is automatically retracted after the puncture and is safely enclosed within the plastic housing.

In particular when dealing with our youngest patients it is imperative that the procedure be as gentle as possible to avoid causing unnecessary anxiety.

The safety mechanism means that the blade/needle is hidden from view, allowing the blood collection to take place in a more relaxed atmosphere for all concerned.



The safety lancets are available with a variety of puncture depths and blade sizes/needle strengths to keep the puncture wound as small as possible.



/ You can find more information about this product in our web-shop.



VACUETTE® QUICKSHIELD

The safety tube holders of the QUICKSHIELD product family are provided with a mechanism that safely keeps the needle in place after activation and thus protects against needlestick injuries.

The VACUETTE® QUICKSHIELD Safety Tube Holder from Greiner Bio-One can optionally be used with VACUETTE® Multiple Use Drawing Needles or with VACUETTE® VISIO PLUS Needles. With our VACUETTE® QUICKSHIELD and VACUETTE® QUICKSHIELD Complete PLUS, the needles are already pre-assembled. The protective shield is at-

tached directly to the holder. The protective shield can be activated with one hand by pressing it towards the needle using a stable surface, alternatively with the finger or thumb. The advantages of this product are the simple handling and effective protection against needlestick injuries and infection.

Without changing the collection technique, the VACUETTE® QUICKSHIELD Safety Tube Holder offers protection against needlestick injuries. The ergonomic design of the holder offers the fingers a better grip during puncture.

- / No change in the blood collection technique
- / More flexibility in terms of user and patient specific requirements due to the rotating protective shield
- / Activation of the protective shield can be felt as well as heard by an audible "click"



VACUETTE® QUICKSHIELD

Safety Tube Holder

The Safety Tube Holder with rotating safety shield is suitable for routine blood collection due to its particularly easy handling. The tube holder can be used with either VACUETTE® Multiple Use Drawing Needles or VACUETTE® VISIO PLUS Needles.

Item no.	Qty. inner / outer
450230	50 / 1,000



VACUETTE® QUICKSHIELD

Complete

VACUETTE® QUICKSHIELD is available as a combination product with a pre-assembled multiple use drawing needle. The individually packaged product is ready for immediate use without any further manual assembly steps.

Item no.	Color code	Needle size	Needle length	Needle length	Qty. inner / outer
450237V1	● green	21 G	1.5 "	38 mm	30 / 360



VACUETTE® QUICKSHIELD

Complete PLUS

VACUETTE® QUICKSHIELD Complete PLUS is available as a combination product with a pre-assembled VISIO PLUS Needle. The transparent viewing window enables visual puncture control. The product is ready for immediate use without any further processing steps.

Item no.	Color code	Needle size	Needle length	Needle length	Qty. inner / outer
450235V1	● green	21 G	1 "	25 mm	30 / 360
450239V1	● green	21 G	1.5 "	38 mm	30 / 360
450226V1	● black	22 G	1 "	25 mm	30 / 360
450228V1	● black	22 G	1.5 "	38 mm	30 / 360



- / You can find more information about this product in our web-shop.



VACUETTE® EVOPROTECT

EVOPROTECT is the next evolution of safety in blood collection.


The semi-automatic release safety mechanism protects the user from the risk of needlestick injuries and makes the daily task of blood collection easier. Because we believe that every needlestick injury is one too many.

The VACUETTE® EVOPROTECT-Safety Blood Collection Set was manufactured with the requirements of the European Directive 2010/32/EU - on preventing sharps injuries in the hospital and healthcare sector - in mind, as well as the OSHA Bloodborne Pathogens Standard. This means that this safety device provides substantial support for safe blood collection or in-

fusion with the same product. User-friendly with an optimized design, the VACUETTE® EVOPROTECT is ergonomically designed and intuitive to use. It has a winged needle, specially developed for one-handed use. The potential needlestick injuries are minimized by the activation of the safety mechanism while the needle is still in the vein.

The triggers of the safety mechanism are located on the side of the product, facilitating one-handed activation. The needle is then automatically retracted to become safely and irreversibly enclosed in the safety shield. A clearly audible click confirms that the safety mechanism is successfully activated. The same device can be used either for taking blood or for a short-term infusion (up to 4 hours).

- / The design and materials of the product are intended for best possible comfort and safety for both user and patient as well as for intuitive application.
- / The transparent view window between the needle and the tubing gives visual control of successful venipuncture
- / The tubing is particularly flexible

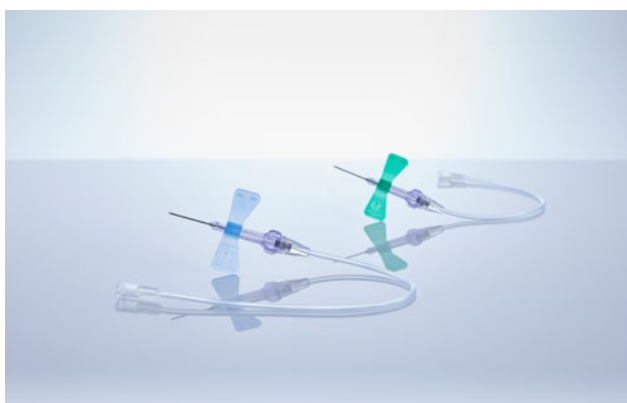


VACUETTE® EVOPROTECT

SAFETY Blood Collection/Infusion Set

This version of the VACUETTE® EVOPROTECT Safety Blood Collection Set is specifically designed for infusion, but may alternatively be used for blood collection. The extra thin wall thicknesses and thus increased flow rates of the 21G and 23G needles enable rapid infusion.

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450121	● blue	23 G	19 mm	7.5 "	19 cm	60 / 1,080

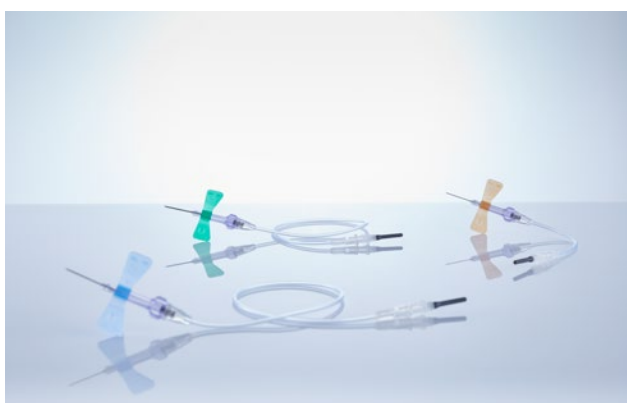

VACUETTE® EVOPROTECT

SAFETY Blood Collection Set + Luer Adapter

VACUETTE® EVOPROTECT is ergonomically designed and intuitive to use. The Luer connection allows the safety product to be used with a holder such as our standard holder or blood culture holder.

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450127	● green	21 G	19 mm	7.5 "	19 cm	60 / 1,080
450128	● green	21 G	19 mm	12 "	30 cm	60 / 1,080
450129	● blue	23 G	19 mm	7.5 "	19 cm	60 / 1,080
450130	● blue	23 G	19 mm	12 "	30 cm	60 / 1,080
450132	● orange	25 G	19 mm	12 "	30 cm	60 / 1,080



VACUETTE® EVOPROTECT

SAFETY Blood Collection Set + Holder

VACUETTE® EVOPROTECT is ergonomically designed and intuitive to use. The safety product is also available with a pre-assembled tube holder, so that you can take the blood sample without any additional previous handling.

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450122	● green	21 G	19 mm	7.5 "	19 cm	30 / 540
450123	● green	21 G	19 mm	12 "	30 cm	30 / 540
450125	● blue	23 G	19 mm	7.5 "	19 cm	30 / 540
450126	● blue	23 G	19 mm	12 "	30 cm	30 / 540



- / You can find more information about this product in our web-shop.



SAFETY BLOOD COLLECTION SET


The SAFETY Blood Collection Set was developed especially for patients with difficult vein conditions, but is equally suitable for any regular blood taking situation, due to the gentle comfort it provides to a patient.

As the protection mechanism is activated whilst the needle is still in the vein, a high level of safety is ensured. The correct activation by the user is indicated by a clearly audible click. The safety blood collection set thus provides the best possible protection against needlestick in-

juries. The visual control of the puncture through a transparent viewing window also increases the puncture safety. The product can be used for blood collection or for infusions for up to five hours.

The SAFETY Blood Collection Set is a sterile single-use winged needle set connected to a flexible tube, with or without a Luer adapter or Luer adapter with holder. The product is ergonomically and individually packed, making it easy to open and minimizing waste.

- / The transparent view window between the needle and the tubing gives visual control of successful venipuncture.
- / Choice of tubing lengths for different applications, featuring minimal memory effect for patient comfort.
- / In-vein activation of safety mechanism to reduce chance of needlestick injuries.

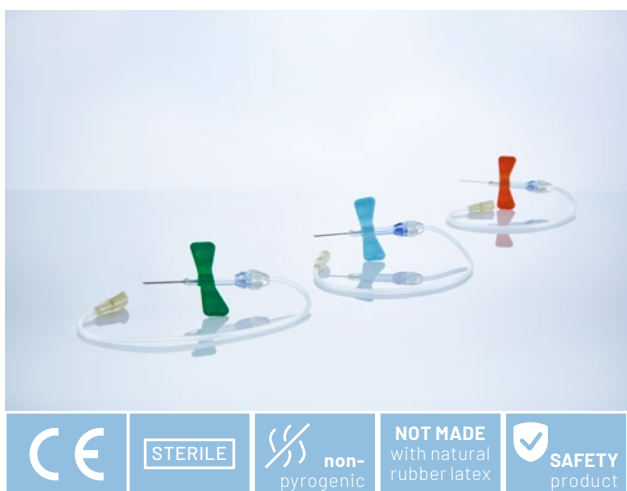

SAFETY Blood Collection/Infusion Set

This version of the product can be used for blood collection by assembling it with, for example, a HOLDEX® or VACUETTE® SAFELINK, or for short-term infusion for up to 5 hours instead.

Safety mechanism: manual activation

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450091V1	● green	21 G	19 mm	7.5 "	19 cm	90 / 1,080
450191V1	● green	21 G	19 mm	12 "	30 cm	90 / 1,080
450092V1	● blue	23 G	19 mm	7.5 "	19 cm	90 / 1,080
450090V1	● orange	25 G	19 mm	7.5 "	19 cm	90 / 1,080



SAFETY Blood Collection Set + Luer Adapter

This version of the SAFETY Blood Collection set has a Luer adapter and can be used with a holder such as our standard tube holder or blood culture holder.

Safety mechanism: manual activation

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450083V1	● green	21 G	19 mm	4 "	10 cm	90 / 1,080
450081V1	● green	21 G	19 mm	7.5 "	19 cm	90 / 1,080
450095V1	● green	21 G	19 mm	12 "	30 cm	90 / 1,080
450084V1	● blue	23 G	19 mm	4 "	10 cm	90 / 1,080
450082V1	● blue	23 G	19 mm	7.5 "	19 cm	90 / 1,080
450096V1	● blue	23 G	19 mm	12 "	30 cm	90 / 1,080
450099V1	● orange	25 G	19 mm	12 "	30 cm	90 / 1,080



SAFETY Blood Collection Set + Holder

This version of the SAFETY Blood Collection set is manufactured with a pre-assembled tube holder and therefore allows an immediate start of the blood collection without assembling with further components.

Safety mechanism: manual activation

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450085V1	● green	21 G	19 mm	7.5 "	19 cm	35 / 420
450160V1	● green	21 G	19 mm	12 "	30 cm	35 / 420
450086V1	● blue	23 G	19 mm	7.5 "	19 cm	35 / 420
450161V1	● blue	23 G	19 mm	12 "	30 cm	35 / 420



SAFETY Blood Collection Set + Blood Culture Holder

This version of the SAFETY Blood Collection set is manufactured with a pre-assembled blood culture holder. The product is suitable for blood collection with common blood culture bottles and subsequently with VACUETTE® Blood Collection Tubes if required.

Safety mechanism: manual activation

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450182	● green	21 G	19 mm	7.5 "	19 cm	24 / 240
450184	● green	21 G	19 mm	12 "	30 cm	24 / 240
450183	● blue	23 G	19 mm	7.5 "	19 cm	24 / 240
450185	● blue	23 G	19 mm	12 "	30 cm	24 / 240



- / You can find more information about this product in our web-shop.



MINICOLLECT® SAFETY LANCETS

Gentle and careful for young & elderly patients.

The pressure activated MiniCollect® Safety Lancets as well as the contact activated MiniCollect® Lancelino Safety Lancets offer a wide range of different puncture depths and thicknesses of blades / needles to keep the puncture wound as small as possible. Due to the safety mechanism, the blade or

needle is automatically retracted after each puncture and safely enclosed in the plastic housing.

- / Safe and hygienic capillary blood collection
- / More precision for easy sample collection
- / Application for glucose tests



MiniCollect®

Lancelino Safety Lancets with Blade

The MiniCollect® Lancelino Safety Lancet with blade is activated simply by pressing the puncture site. Due to the round contact surface a precise placement is possible. The product with blade is ideal for the collection of larger blood volumes.

Safety mechanism: automatic activation

Item no.	Color code	Blade size	Penetration depth	Qty. inner / outer
450524	● yellow	0.8 mm	2 mm	200 / 2,000



MiniCollect®

Lancelino Safety Lancets with Needle

The MiniCollect® Lancelino Safety Lancet with needle is activated simply by pressing the puncture site. Due to the round contact surface a precise placement is possible. The product with needle is perfectly suited for the collection of smaller capillary blood volumes.

Safety mechanism: automatic activation

Item no.	Color code	Needle size	Penetration depth	Qty. inner / outer
450521	● lavender	25 G	1.5 mm	200 / 2,000
450522	● blue	21 G	1.8 mm	200 / 2,000
450523	● green	21 G	2.4 mm	200 / 2,000



MiniCollect®

Safety Lancets with Blade

MiniCollect® Safety Lancets enable safe and hygienic capillary blood collection. The blade only releases when the release button is pressed. The product with blade is suitable for the collection of small to large capillary blood volumes.

Safety mechanism: automatic activation

Item no.	Color code	Blade size	Penetration depth	Qty. inner / outer
450427	● pink	1.5 mm	1 mm	200 / 2,000
450428	● green	1.5 mm	1.5 mm	200 / 2,000
450429	● blue	1.5 mm	2 mm	200 / 2,000




MiniCollect®

Safety Lancets with Needle

MiniCollect® Safety Lancets enable safe and hygienic capillary blood collection. The needle only releases when the release button is pressed. The product with needle is best suited for the collection of single blood drops.

Safety mechanism: automatic activation

Item no.	Color code	Needle size	Penetration depth	Qty. inner / outer
450438	● lavender	28 G	1.25 mm	200 / 2,000
450439	● orange	23 G	2.25 mm	200 / 2,000

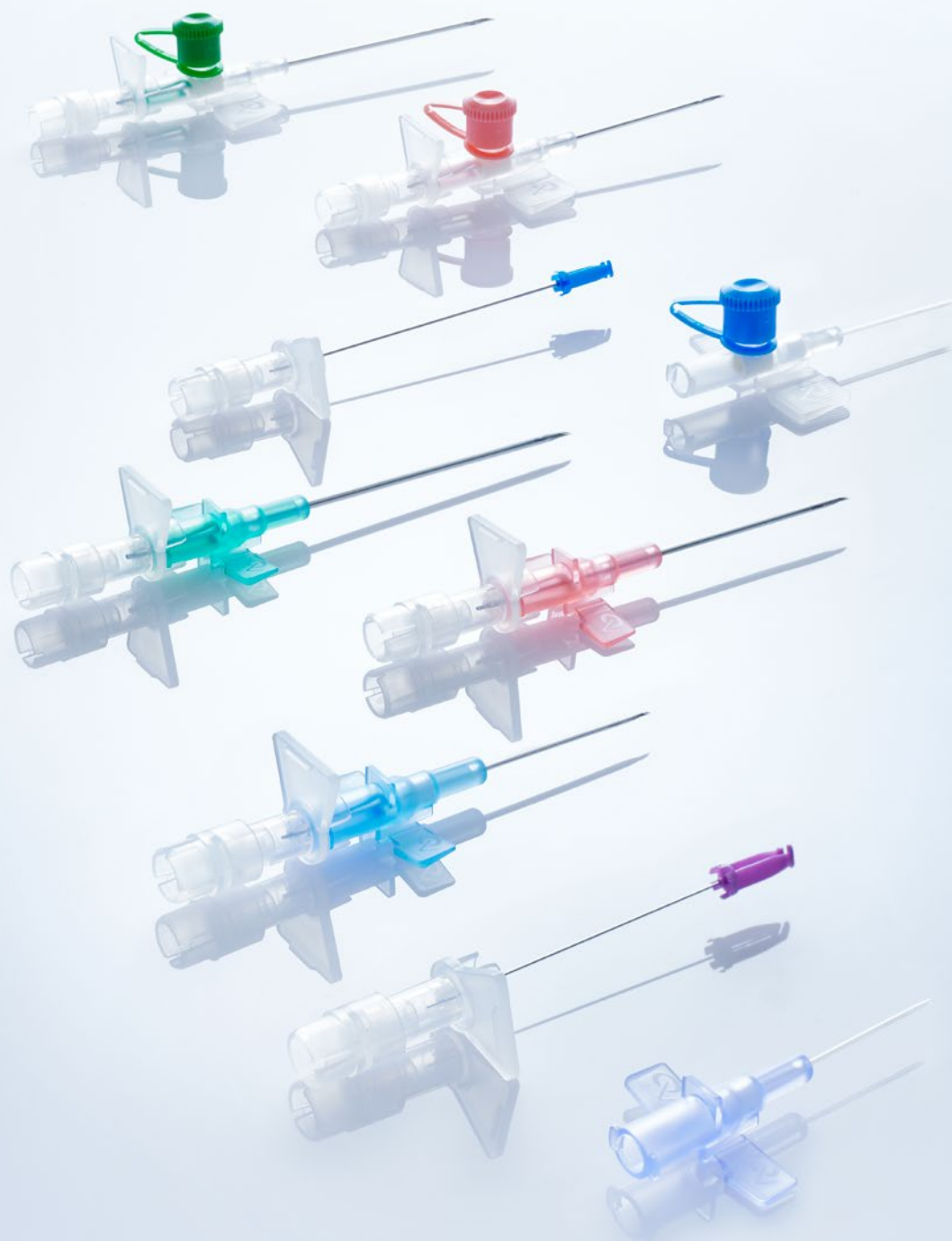

Our IV catheters have excellent insertion properties and are user-friendly, safe, sharp and with a consistent quality, so you can focus less on worrying about needlestick injuries and more on giving the best care to your patients.

IV CATHETERS

/	CLiP® Automatic Safety IV Catheter	102
	CLiP® Neo Safety IV Catheter	103
	CLiP® Winged Safety IV Catheter.....	103
	CLiP® Ported Automatic Safety IV Catheter...	104



/ You can also find out more about our safety catheters at www.gbo.com.



CLiP® AUTOMATIC SAFETY IV CATHETER

CLiP Automatic Safety IV Catheter enables a safe and successful insertion.

Our IV catheters are designed with the healthcare professional's safety and working environment in mind. Our focus is to provide products that are user-friendly, safe, sharp and with a consistent quality. CLiP® safety IV catheters are

available in many variations to meet demanding physical conditions and different clinical needs. CLiP® Neo is ideal for the smallest patients - newborns and infants - or elderly people. It was developed especially for small and fragile veins. CLiP® Winged and CLiP® Ported are available in seven gauge sizes, two different material and many different lengths.

- / Automatic safety
- / Excellent insertion properties
- / Comprehensive product portfolio
- / Slim packaging design for streamlined storage space



CLiP® Neo Safety IV Catheter

CLiP® Neo is an automatic safety IV catheter designed for small and fragile veins found in the neonatal, pediatric or elderly patient groups. The product is available in 24 and 26 gauges to cater for demanding physical conditions. The optimized needle notch provides instant blood response.

Item no.	Catheter material	Needle size	Catheter dimension	Color code	Catheter flowrate	Qty. inner / outer
NW241901	FEP	24 G	0.7 x 19	yellow	22	500 / 1,000
NW261901	FEP	26 G	0.6 x 19	violet	15	500 / 1,000



CLiP® Winged Safety IV Catheter


CLiP® Winged is an automatic safety IV catheter available in a comprehensive range (14-24 gauge). CLiP® Winged is suitable for high pressure injections (up to 305 psi) and has wing holder to enable multiple gripping alternatives.

Item no.	Catheter material	Needle size	Catheter dimension	Color code	Catheter flowrate	Qty. inner / outer
VW144511	PUR	14 G	2.0 x 45	orange	290	500 / 1,000
VW164511	PUR	16 G	1.7 x 45	gray	200	500 / 1,000
VW174511	PUR	17 G	1.5 x 45	white	140	500 / 1,000
VW183211	PUR	18 G	1.2 x 32	green	110	500 / 1,000
VW184511	PUR	18 G	1.2 x 45	green	100	500 / 1,000
VW203211	PUR	20 G	1.0 x 32	pink	64	500 / 1,000
VW222511	PUR	22 G	0.8 x 25	blue	38	500 / 1,000
VW241911	PUR	24 G	0.7 x 19	yellow	22	500 / 1,000
VW144501	FEP	14 G	2.0 x 45	orange	290	500 / 1,000

IV catheters

CLiP® Automatic Safety IV Catheter

Item no.	Catheter material	Needle size	Catheter dimension	Color code	Catheter flowrate	Qty. inner / outer
VW164501	FEP	16 G	1.7 x 45	● gray	200	500 / 1,000
VW174501	FEP	17 G	1.5 x 45	○ white	140	500 / 1,000
VW183201	FEP	18 G	1.2 x 32	● green	110	500 / 1,000
VW184501	FEP	18 G	1.2 x 45	● green	100	500 / 1,000
VW203201	FEP	20 G	1.0 x 32	● pink	64	500 / 1,000
VW222501	FEP	22 G	0.8 x 25	● blue	38	500 / 1,000

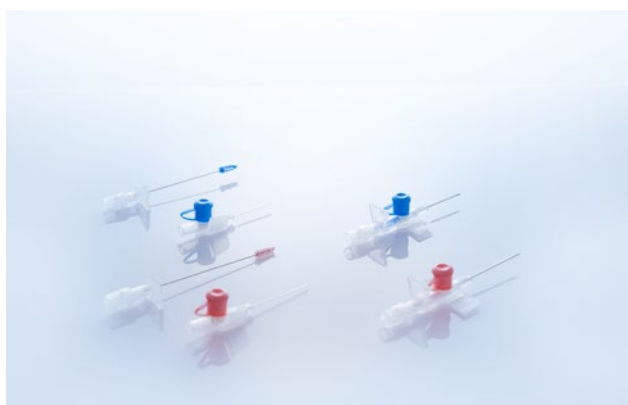

CLiP® Ported Automatic Safety IV Catheter

CLiP® Ported is an automatic safety IV catheter with injection port that is suitable for high pressure injections (up to 305 psi). The product is available in a comprehensive range to meet demanding physical conditions and different clinical needs.



Item no.	Catheter material	Needle size	Catheter dimension	Color code	Catheter flowrate	Qty. inner / outer
VP144511	PUR	14 G	2.0 x 45	● orange	290	500 / 1,000
VP164511	PUR	16 G	1.7 x 45	● gray	200	500 / 1,000
VP174511	PUR	17 G	1.5 x 45	○ white	140	500 / 1,000
VP183211	PUR	18 G	1.2 x 32	● green	110	500 / 1,000
VP184511	PUR	18 G	1.2 x 45	● green	100	500 / 1,000
VP203211	PUR	20 G	1.0 x 32	● pink	64	500 / 1,000
VP202511	PUR	20 G	1.0 x 25	● pink	68	500 / 1,000
VP222511	PUR	22 G	0.8 x 25	● blue	38	500 / 1,000
VP241911	PUR	24 G	0.7 x 19	● yellow	22	500 / 1,000
VP144501	FEP	14 G	2.0 x 45	● orange	290	500 / 1,000
VP164501	FEP	16 G	1.7 x 45	● gray	200	500 / 1,000
VP174501	FEP	17 G	1.5 x 45	○ white	140	500 / 1,000
VP183201	FEP	18 G	1.2 x 32	● green	110	500 / 1,000
VP184501	FEP	18 G	1.2 x 45	● green	100	500 / 1,000
VP202501	FEP	20 G	1.0 x 25	● pink	68	500 / 1,000
VP203201	FEP	20 G	1.0 x 32	● pink	64	500 / 1,000
VP222501	FEP	22 G	0.8 x 25	● blue	38	500 / 1,000
VP241901	FEP	24 G	0.7 x 19	● yellow	22	500 / 1,000



Greiner Bio-One offers a variety of blood collection needles and blood collection sets for your individual needs. The products are mainly used for routine blood collection and are intended for single use only.

SHARPS

/ Blood Collection Needles	108
VACUETTE® VISIO PLUS Needles.....	109
VACUETTE® Multiple Use Drawing Needles....	109
/ Blood Collection Sets	110
Blood Collection Set + Luer Adapter	111
Blood Collection Set + Holder	111



- / We also offer blood collection needles for veterinary use. Please refer to our website.
- / You can find more information about this product in our web-shop.



BLOOD COLLECTION NEEDLES

At Greiner Bio-One you will find a wide range of blood collection needles for human use but also for veterinary use.

The range of product options ensure that the right products for your daily routine works are available. Blood collection needles are intended for single use only.

Our range includes the Multiple Use Drawing Needles as well as the VACUETTE® VISIO PLUS Needles with optical puncture control. A variety of options are available in terms of dimension and needle length. For veterinary use, Single Sample Needles are also on offer.

- / Sterile
- / Color-coded needle dimensions
- / Available in different lengths
- / VISIO PLUS Needle with optical puncture control



VACUETTE® VISIO PLUS Needles

The VACUETTE® VISIO PLUS Needle with optical venipuncture control can simplify blood collection routine for phlebotomists.

Item no.	Color code	Needle size	Needle length	Needle length	Qty. inner / outer
450042	● green	21 G	1 "	25 mm	2,000 / 6,000
450043	● black	22 G	1 "	25 mm	2,000 / 6,000
450040	● green	21 G	1.5 "	38 mm	2,000 / 6,000
450041	● black	22 G	1.5 "	38 mm	2,000 / 6,000



VACUETTE® Multiple Use Drawing Needles

Blood Collection Needles are predominantly used in routine blood collection and are exclusively single-use. The color of the needle protective cap gives a visual indication of the needle gauge size.

Item no.	Color code	Needle size	Needle length	Needle length	Qty. inner / outer
450071	● black	22 G	1 "	25 mm	2,000 / 6,000
450072	● green	21 G	1 "	25 mm	2,000 / 6,000
450073	● yellow	20 G	1 "	25 mm	2,000 / 6,000
450075	● black	22 G	1.5 "	38 mm	2,000 / 6,000
450076	● green	21 G	1.5 "	38 mm	2,000 / 6,000
450077	● yellow	20 G	1.5 "	38 mm	2,000 / 6,000



/ You can find more information about this product in our web-shop.



BLOOD COLLECTION SETS

Blood collection sets, also known as winged needles or butterfly, are a popular product for blood collection and are available in various types.

These products are primarily for use on patients with difficult vein conditions. Blood Collection Sets are available with dif-

ferent tubing lengths and needle gauges.

The sets are available pre-assembled with Luer Adapter or Luer Adapter and Holder.

Blood collection sets are intended to be used for blood collection.

/ **Blood Collection Set + Luer Adapter**

/ **Blood Collection Set + Holder**



Blood Collection Set + Luer Adapter

The Blood Collection Set + Luer Adapter is a sterile single-use winged needle connected to a flexible tubing. This product version can be used with a holder such as our standard tube holder or blood culture holder.

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450064	● green	21 G	19 mm	4 "	10 cm	50 / 1,000
450065	● green	21 G	19 mm	12 "	30 cm	50 / 1,000
450066	● green	21 G	19 mm	7.5 "	19 cm	50 / 1,000
450089	● blue	23 G	19 mm	7.5 "	19 cm	50 / 1,000
450159	● blue	23 G	19 mm	12 "	30 cm	50 / 1,000



Blood Collection Set + Holder

The Blood Collection Set + Holder is a sterile single-use winged needle connected to a flexible tubing. This version is manufactured with a pre-assembled holder, enabling an immediate start of the blood collection without assembling of further components.

Item no.	Color code	Needle size	Needle length	Tubing length	Tubing length	Qty. inner / outer
450373	● blue	23 G	19 mm	7.5 "	19 cm	24 / 240



During the development of all accessories, special attention is paid to ensuring that the products absolutely meet the requirements of the end users. The product range includes various tube holders, Luer adapters, tourniquets, ID snap rings and much more.

BLOOD COLLECTION ACCESSORIES

/ Tube Holder	114
Standard Tube Holder	115
HOLDEX® Single-Use Holders	115
VACUETTE® SAFELINK	116
VACUETTE® Blood Transfer Unit	116
VACUETTE® Blood Culture Holder	117
 / VACUETTE® SECONDARY Tubes.....	118
VACUETTE® SECONDARY Tubes SIMPLEX	119
VACUETTE® SECONDARY Tubes MULTIPLEX ..	119
 / Luer Adapter	120
VACUETTE® Luer Adapter	121
 / VACUETTE® Tourniquets.....	122
VACUETTE® Tourniquets	123
VACUETTE® Disposable Tourniquets	123
VACUETTE® EasyClean Tourniquets	124
 / VACUETTE® ID Snap Rings	125
VACUETTE® ID Snap Rings	126
 / Blood Collection Tray	127
Blood Collection Tray	128



/ You can find more information about this product in our web-shop.



TUBE HOLDER

Simple handling and safety are of the highest importance during daily routine procedures.

In the development of tube holders, special attention is paid to ensuring that the products meet the requirements of the users. To ensure optimal performance,

Tube holders should be used exclusively with blood collection accessories from Greiner Bio-One, as the components are designed for optimal compatibility with one another.

- / Standard Tube Holder
- / HOLDEX® Single-Use Holders
- / SAFELINK
- / Blood Transfer Unit
- / Blood Culture Holder



Standard Tube Holder

The standard tube holder should be used in the daily blood collection routine. It should be used exclusively with VACUETTE® Multiple Use Drawing Needles or VACUETTE® VISIO PLUS Needles, as both components are designed for optimal compatibility with one another.

Item no.	Qty. inner / outer
450209	100 / 2,000



STERILE

HOLDEX® Single-Use Holders

The HOLDEX® single-use plastic holder features an off-center luer connector. This special design of the holder provides a flatter angle for venipuncture. The luer connector of the HOLDEX® single-use holders allows quick and easy connection with blood collection accessories with female luer connector, whilst the practical plastic needle of the Plastic Cannula HOLDEX® enables safe sample collection from needleless ports with rubber membrane.

Item no.	Description	Qty. inner / outer
450241	HOLDEX® Single-Use Holder PP, single-packed, not made with natural rubber latex	100 / 800
450263	HOLDEX® Single-Use Holder PP, single-packed	100 / 800
450216	Plastic Cannula HOLDEX®, single-packed	1 / 1,000



STERILE

NOT MADE
with natural
rubber latex

VACUETTE® SAFELINK

The innovative product allows easy and quick blood sampling when used with medical devices that have a standardized female luer lock connection. VACUETTE® SAFELINK can be combined with all standardized female luer lock connections for blood sampling.

Description: Holder with male luer lock, single-packed

Item no.	Qty. inner / outer
450210	100 / 800



IVD

STERILE

VACUETTE® Blood Transfer Unit

The blood transfer unit (BTU) is for the safe, needle-less transfer of samples from a syringe into an evacuated blood collection tube. The female luer connection allows for safe sample transfer from standardized syringes.

Description: PP, single-packed

Item no.	Qty. inner / outer
450225	100 / 800



NOT MADE
with natural
rubber latex

VACUETTE® Blood Culture Holder

The blood culture holder is a non-sterile single-use device designed for filling commonly used blood culture bottles, in addition to VACUETTE® Blood Collection Tubes.

Description: VACUETTE® Blood Culture Holder, non-sterile

Item no.	Qty. inner / outer
450181	40 / 320



- / VACUETTE® Secondary Tubes can be easily closed with our Snap Caps.
- / You can find more information about this product in our web-shop.



VACUETTE® SECONDARY TUBES

VACUETTE® SECONDARY Tubes support laboratory automation with their versatile application.

Laboratory automation continues to progress and consequently sample throughput is on the rise. In order to encourage efficient procedures in the laboratory, Greiner Bio-One has developed the VACUETTE® SECONDARY Tubes. The greatest advantage of this product line is its versatility.

The non-ridged version of the

SIMPLEX tubes is already packed in alignment so they can be transferred directly to the storage container (plastic or metal box) of the pre-analytical sample processing system. SIMPLEX tubes are manufactured from PET.

MULTIPLEX tubes are available in sealed plastic bags for manual application. The PREMIUM tubes are non-sterile and are available in PET or PP, with or without caps.

- / Aliquoting
- / Pooling
- / Storage and/or transport
- / Centrifugation
- / Freezing or heating



VACUETTE® SECONDARY Tubes SIMPLEX

The non-ridged version of the SIMPLEX tubes is already packed in alignment so they can be transferred directly to the storage container (plastic or metal box) of the pre-analytical sample processing system. SIMPLEX tubes are manufactured from PET.

Cap design: pull cap

Item no.	Description	Qty. inner / outer
459020	RD PET, 13x75 mm without cap, non-ridged	496 / 1,984
459021	BC PET, 13x75 mm without cap, non-ridged	496 / 1,984

- / 459020: designed for Roche PVT instrumentation (cobas p612)
- / 459021: designed for Beckman Coulter (Olympus) instrumentation (AutoMates 1250 + 2550)
- / VACUETTE® SECONDARY Tubes SIMPLEX are non-sterile and can be used at temperatures from -80 °C up to +50 °C.



VACUETTE® SECONDARY Tubes MULTIPLEX

MULTIPLEX tubes are available in sealed plastic bags for manual application. The PREMIUM tubes are non-sterile and are available in PET or PP, with or without caps.

Cap design: safety screw cap

Item no.	Description	Qty. inner / outer
459000	PP, 13x75 mm with white screw closure, PREMIUM	500 / 2,000
459001	PP, 13x75 mm with orange screw closure, PREMIUM	500 / 2,000
459010	PP, 13x75 mm without cap, PREMIUM	500 / 2,000
459011	PET, 13x75 mm without cap, PREMIUM	500 / 2,000

- / The PET tubes can be used at temperatures between -80 °C and +50 °C, the PP tubes from -80 °C up to +121 °C.
- / In combination with screw caps, the MULTIPLEX tubes fulfill the IATA guidelines (95 kPa) for secure sample transport and are compliant with the CLSI Standard ATU001-A.



/ You can find more information about this product in our web-shop.



LUER ADAPTER

The Luer adapter is a standardized one-size-connection system for needle devices.

The Luer adapter is manufactured with a standardized Luer part which allows for quick and easy assembly with the chosen

component device.

Luer adapters have a male Luer and are suitable for use with products like SAFETY Blood Collection Sets and other needle devices with a female Luer connector.

- / Easy to connect
- / Luer Adapter with male Luer
- / Compatible with, for example, SAFETY Blood Collection Sets or Luer needles



CE

STERILE

 non-
pyrogenic
NOT MADE
with natural
rubber latex**VACUETTE® Luer Adapter**

The VACUETTE® Luer Adapter is for use in venous blood collection with attached winged needle or luer needle. The product is a single-use device and is fitted with a safety valve at one end.

Item no.	Color code	Needle size	Qty. inner / outer
450070	● royal blue	20 G	2,000 / 6,000



/ You can find more information about this product in our web-shop.



VACUETTE® TOURNIQUETS

VACUETTE® tourniquets round off the Greiner Bio-One product range for successful blood collection.

Greiner Bio-One has a choice of latex-free single-use or reusable tourniquets. The VACUETTE® Tourniquet can be used several times and features a buckle enabling easy one-handed adjustment and release. Also availa-

ble in a version especially for children with attractive, colorful designs.

The VACUETTE® Disposable Tourniquet is for single-use only and comes perforated on a roll in a practical dispenser box.

There is also a suitable product for your younger patients. The VACUETTE® Tourniquet is available in very special designs as a version for children.

- / Latex-free
- / Easy handling
- / Single-use and multiple-use tourniquets
- / Special designs for children's tourniquets



NOT MADE
with natural
rubber latex

VACUETTE® Tourniquets

The VACUETTE® Tourniquet is designed to help find a venipuncture site and is characterized by one-handed use with a quick release belt. This product is intended for multiple use. Also available in special designs for children.

Length: 500 mm

Item no.	Description	Qty. inner / outer
840050	non-sterile	1 / 100
840051	for children, non-sterile	- /100




NOT MADE
with natural
rubber latex

VACUETTE® Disposable Tourniquets

The disposable tourniquets from Greiner Bio-One are characterized by easy handling. The VACUETTE® disposable tourniquet is latex free, soft on the skin and intended for single-use only.

Description: VACUETTE® Disposable Tourniquet, synthetic rubber, non-sterile, Length: 480 mm

Item no.	Qty. inner / outer
840053	500 / 2,500

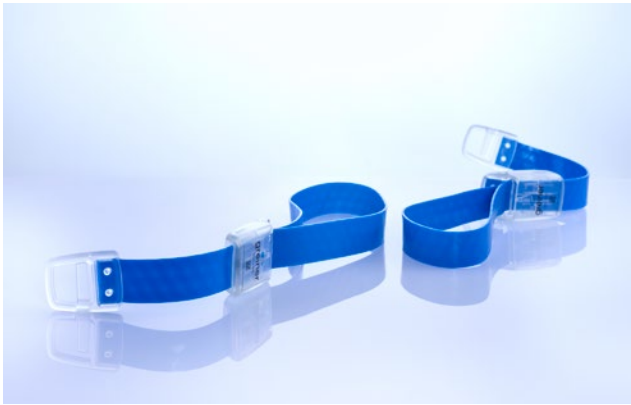

VACUETTE® EasyClean Tourniquets

The VACUETTE® EasyClean Tourniquet is made out of silicone and has a specially treated surface for easy cleaning between uses with simple wipe disinfection. The VACUETTE® EasyClean features the easy application of a classic buckle tourniquet, including one-handed release and adjustment but with added hygienic safety due to the easy wipe disinfection. Tested for use up to 11,200 times.



Item no.	Qty. inner / outer
840052	1 / 80



- / You can find more information about this product in our web-shop.



VACUETTE® ID SNAP RINGS

VACUETTE® ID Snap Rings are a practical accessory product.

With the help of the ID Snap Rings, tubes can be individually differentiated at first glance, for

example for special analysis purposes in standard processes. The VACUETTE® ID Snap Rings fit all VACUETTE® 13 mm tube caps as well as all 16 mm tube caps.

- / Available in various colors
- / For easy, individual identification of 13 mm and 16 mm tubes



VACUETTE® ID Snap Rings

VACUETTE® ID Snap Rings are available in all five colors and fit all VACUETTE® 13 mm and 16 mm tubes. With the Snap Rings, you can easily and quickly mark tubes individually.

Item no.	Product color	Qty. inner / outer
458711	● orange	500 / 4,000
458712	● green	500 / 4,000
458713	● pink	500 / 4,000
458714	● light blue	500 / 4,000
458715	● grey	500 / 4,000



- / You can find more information about this product in our web-shop.



BLOOD COLLECTION TRAY

- / Perfectly organized with the blood collection tray
- / Lightweight
- / Easy to clean
- / Materials and products always at hand
- / Compatible with sharps disposal boxes from 0.6 to 3 liters

With the Compact Blood Collection Tray you are optimally equipped for blood collection.

The tray is used to store all materials and products required for blood collection, so that every-

thing is always at hand. The blood collection tray is lightweight and easy to clean.

The tray is suitable for our sharps disposal boxes from 0.6 to 3 liters (see also chapter «Disposal»).



Blood Collection Tray

With the Blood Collection Tray you are perfectly organized for blood collection and always have your materials and products at hand. The tray is easy to clean, lightweight and compatible with sharps disposal boxes from 0.6 to 3 liters.

Description: white, for 0.6 / 1 / 1.8 / 3 liter Sharps Disposal Containers

Item no.	Qty. inner / outer
450343	1 / 22



Needlestick injuries are a particular risk during the disposal of sharp or pointed objects. To reduce the risk for the user to a minimum, only containers specifically for this purpose should be used. Furthermore, the safety mechanism of the safety product should be activated prior to disposal.

DISPOSAL

/ Disposal Containers 132

 BioSecure Plus Disposal Boxes 133

 Biocompact Disposal Boxes 133

 Biogrip Disposal Boxes 134



/ You can find more information about this product in our web-shop.



DISPOSAL CONTAINERS

Disposal containers are part of our product line and enable the safe disposal of medical consumables. All containers are of penetration-proof plastic and feature a temporary and permanent closure.

Blood collection accessories

as well as all other pointed and sharp objects must be disposed of in sharps-proof boxes. All disposal containers from Greiner Bio-One have a large opening for the disposal of sharp and pointed objects and are available in sizes from 0.6 liters up to 50 liters.

- / Penetration-proof plastic
- / Temporary and irreversible closure
- / Comprehensive product range in different sizes



BioSecure Plus Disposal Boxes

BioSecure Plus containers, available in 2 fill volumes, are manufactured from 40% post-consumer recycled (PCR) content. The wide lid opening enables safe and easy disposal of larger items like blood collection sets. The lid design also includes three sharps removal tabs, including one specifically for insulin needles.

Item no.	Product color	Nominal volume	Qty. inner / outer
450330V1	● yellow	0.6 l	1 / 75
450331V1	● yellow	1 l	1 / 75



Biocompact Disposal Boxes

Biocompact disposal boxes are available in sizes between 1.8 liters and 10 liters. The container is suitable for the disposal of contaminated pointed, sharp and biohazardous materials. The product is puncture resistant.

Item no.	Product color	Nominal volume	Qty. inner / outer
450332	● yellow	1.8 l	1 / 30
450333	● yellow	3 l	1 / 30
450334	● yellow	5 l	1 / 30
450335	● yellow	10 l	1 / 20



Biogrip Disposal Boxes

Biogrip disposal boxes are available in 30 liter and 50 liter sizes. The container is suitable for the disposal of contaminated pointed, sharp and biohazardous materials. The product is puncture resistant.

Item no.	Product color	Nominal volume	Qty. inner / outer
450338	● yellow	30 l	- /1
450339	● yellow	50 l	- /1



Adequate and proper sample transport helps to keep the samples cool and to transport them safely. This is an extremely important influencing factor for getting correct analysis results. The products of the VACUETTE® Transport Line enable the safe transport of biological substances (UN 3373) of the category B – from collection to analysis in the laboratory.

SAMPLE TRANSPORT

/ VACUETTE® Transport Line.....	138
VACUETTE® Transport Boxes	139
VACUETTE® Transport Containers	139
/ VACUETTE® Transport Line	
Accessories.....	140
Transport Bags	141



- / You can also find suitable transport bags under accessories VACUETTE® transport line.
- / You can find more information about this product in our web-shop.



VACUETTE® TRANSPORT LINE

Better safe than sorry. Sample transport has never been easier.

With the VACUETTE® Transport Line, your valuable samples are transported safely and kept cool. The range includes transport boxes and containers for your special sample transport requirements.

The combination of primary, secondary and outer packaging (3rd packaging) from our VACUETTE® transport line corresponds to ADR/RID (packa-

ging instructions P650), IATA-DG, ICAO-TI and IMDG-Code.

It's not all just about regulatory compliance and sample quality. The design of our transport solutions also ensures workflow efficiency and user comfort. Quick and easy loading and unloading of samples, as well as the extreme robustness all contribute to reducing the challenges of everyday routine for healthcare professionals.

- / Safe sample transport
- / Tested and certified by an independent, accredited institute
- / Large variety of accessories



VACUETTE® Transport Boxes

The VACUETTE® transport boxes made of polypropylene (PP) are stackable and suitable for transporting up to 40 tubes. The box is liquid-tight and has a click closure for fast and convenient opening and closing. The foam insert is suitable for tube diameters of 9-16 mm.

Item no.	Description	Qty. inner / outer
472001	VACUETTE® Transport Box (VTB), incl. Foam Insert for 40 Tubes, with Transport Carton HK0190	1 / 10
472040	VACUETTE® Transport Box (VTB), incl. Foam Insert for 40 Tubes	1 / 10
472010	Foam Insert for VACUETTE® Transport Box (VTB), for 40 Tubes (all sizes)	- / 168
472016	Cooler Pack for VACUETTE® Transport Box (VTB)	- / 72
472011	Rack for Blood Collection Tubes, 10 positions, autoclavable	1 / 100
472015	Absorbent Pad for VTB / VTC	- / 2,500



VACUETTE® Transport Containers

Our VACUETTE® Transport Container is a handy container made of polypropylene (PP) for transporting up to 12 tubes. The container is closed with a screw cap and is liquid-tight. The foam insert is suitable for tube diameters of 9-16 mm.

Item no.	Description	Height	Width	Qty. inner / outer
800105	VACUETTE® Transport Container (VTC), incl. Foam Insert for 12 Tubes, with Transport Carton HK0259	128.7 mm	113 mm	1 / 12
800110	VACUETTE® Transport Container (VTC), incl. Foam Insert for 12 Tubes, without Transport Carton HK0259	128.7 mm	113 mm	1 / 12
800101	Foam insert for VACUETTE® Transport Container (VTC), for 12 tubes	70 mm	107 mm	- / 336
800102	Container volume 60 ml with absorbent pad and screw cap, without symbol	127 mm	28 mm	- / 300
472017	Cooler Pack for VACUETTE® Transport Container (VTC)	12 mm	107 mm	- / 150



- / Test reports can be downloaded from the download center on our website: www.gbo.com
- / You can find more information about this product in our web-shop.



VACUETTE®

TRANSPORT LINE ACCESSORIES

Our extensive range of accessory products optimally rounds off the VACUETTE® transport line.

In addition to our transport boxes and our transport container, we also provide a wide range of isothermal bags as an alternative to transport cartons, in which you can carry transport boxes or containers.

The VACUETTE® transport line also includes cooling packs, cooling pads, one foam insert each for the transport box and the transport container and inserts as well as labels for the outer packaging.

Greiner Bio-One offers you all you need for safe sample transport.

- / Safe sample transport
- / Large variety of accessories
- / Cooling options
- / Racks, foam insert, absorbent inlays



Transport Bags

The Isotherm transport bags are available in different sizes for the transport of several boxes or containers. The transport bags are robust, light weight and help to keep the temperature low due to the thermally insulating outer material.

Item no.	Description	Height	Length	Width	Product color	Qty. inner / outer
472020	for 1 piece VACUETTE® Transport Box (VTB), small	250 mm	300 mm	170 mm	● royal blue	1 / 16
472030	for 2 pieces VACUETTE® Transport Boxes (VTB), medium	250 mm	320 mm	300 mm	● royal blue	1 / 8
472023	for 3 pieces VACUETTE® Transport Boxes (VTB), medium	250 mm	480 mm	300 mm	● royal blue	1 / 4
472024	for 4 pieces VACUETTE® Transport Boxes (VTB), big	250 mm	640 mm	300 mm	● royal blue	1 / 4
472022	for 1 piece VACUETTE® Transport Container (VTC), small	160 mm	145 mm		● royal blue	1 / 50



At Greiner Bio-One, you will find numerous products and accessories for general laboratory requirements in top quality. Our product range includes Snap Caps, carrier tubes and caps for MiniCollect® and much more.

GENERAL LABORATORY EQUIPMENT

/ Snap Caps	144
Snap Caps	145
/ MiniCollect® Accessories Tubes	146
MiniCollect® Carrier Tube	147
MiniCollect® Extra Cap.....	147



/ You can find more information about this product in our web-shop.



SNAP CAPS

Snap Caps from Greiner Bio-One are mainly used to: recap open sample collection tubes.

The Snap Caps are caps, available in a variety of different colors

and can be used to recap 13mm non-ridged and PREMIUM PET tubes, both VACUETTE® and Simplex BC. Not suitable for transport.

- / Available in various colors
- / For VACUETTE® Tubes with 13mm diameter
- / Made of polyethylene (PE)
- / Not suitable for transport



Snap Caps

Snap Caps from Greiner Bio-One are available in many different colors and are used to recap open sample collection tubes. The product is suitable for VACUETTE® tubes with a diameter of 13 mm.

Cap design: pull cap

Item no.	Product color	Qty. inner / outer
371516	○ white	1,000 / 5,000
371517	● green	1,000 / 5,000
371518	● yellow	1,000 / 5,000
371519	● royal blue	1,000 / 5,000
371522	● red	1,000 / 5,000
371523	● lavender	1,000 / 5,000
371527	● pink	1,000 / 5,000



- / You can find more information about this product in our web-shop.



MINICOLLECT® ACCESSORIES TUBES

- / Available as transparent or amber (light protection) carrier tube
- / Processing and analysis of capillary and venous blood samples in the same tube dimension
- / PREMIUM Tubes
- / Extra caps for MiniCollect® tubes

MiniCollect® is the ideal system for the collection, transport and processing of capillary blood. When our carrier tube meets MiniCollect®, it is the perfect combination.

MiniCollect® tubes can be threaded into a PREMIUM carrier tube. This combination corresponds to the dimensions of a standard 13 x 75mm tube.

This means that it is not necessary to change the settings on the analyzer for many device types. Capillary and venous blood samples can be processed and analyzed in the same tube dimension.

In addition to our MiniCollect® carrier tubes, the Greiner Bio-One range also includes extra caps for the MiniCollect® tubes.



MiniCollect® Carrier Tube

MiniCollect® tubes can be threaded into a PREMIUM carrier tube. This combination corresponds to the dimensions of a standard 13 x 75mm tube. The carrier tubes are available in transparent and amber.

Item no.	Product color	Thread type	Qty. inner / outer
450417	○ clear	PREMIUM	100 / 500
450418	● amber	PREMIUM	100 / 500



MiniCollect® Extra Cap

The accessories for MiniCollect® tubes include spare caps to enable quick and easy recapping of MiniCollect® blood collection tubes after analysis.

Cap design: pull cap

Item no.	Product color	Qty. inner / outer
370267B	○ white	- / 9,008



VeinViewer® packs the latest advances in superior vascular imaging into a mobile package. VeinViewer® can help deliver faster venipuncture procedures on your difficult patients through greatly improved first stick success. It gives you revolutionary HD views in realtime, so you have the best view of vasculature procedures pre-, during and post-access. MaxDepth allows a view of peripheral veins up to 10mm deep and blood patterns up to 15mm beneath the surface of the skin and CathCompass™ is an integrated reference graphic to help choose the right IV catheter size.

DEVICES

/ VeinViewer® 150

 VeinViewer® Flex 151

 VeinViewer® Vision2 151



/ You can find more information about this product in our web-shop.



VEINVIEWER®

With the near infrared light technology of VeinViewer® you benefit before, during and after venipuncture. The device makes it possible to see through the patient's skin and identify vessels. The CathCompass™ feature provides an integrated reference graphic for IV catheter sizes projected onto the skin.

The VeinViewer® uses harmless near-infrared (NIR) light which is directed towards the patient's skin. Hemoglobin in the blood absorbs the NIR and the surrounding tissue reflects it back to the

VeinViewer® device, where the data is processed into an image, color is added and the image sent back to the skin to provide a real time visualization of the blood vessels and patterns up to 10mm below the skin.

HD technology with digital full field help provide the clinician with clear, sharp imaging to support routine healthcare work. VeinViewer® is available as a portable device VeinViewer® Flex or as a mobile device with an arm with a range of up to 1.35 meter as VeinViewer® Vision2.

- / Improved first stick success
- / Improved patient satisfaction
- / Improved safety



VeinViewer® Flex

VeinViewer® Flex is a hand-carried device which is also suitable for use in emergency situations. The product can be made static by using the flexible arm provided for use anywhere, especially if space is limited. It is easy to use and impresses with quality, durability and maximum mobility.

Description: Contains 2 pcs. Lithium Ion batteries

Item no.	Qty. inner / outer
890001	- /1

/ Only available on selected markets



VeinViewer® Vision2

VeinViewer® Vision2 is a mobile imaging unit with an arm reaching up to 1.35 meter. It is attached to a robust trolley for easy transportation between beds or phlebotomy chairs.

Description: Contains 1 pc. Lithium Ion battery

Item no.	Qty. inner / outer
890004	- /1

/ Only available on selected markets



Urine, as a metabolic byproduct, can tell a lot about a person's health status and any changes observed may indicate a certain medical condition. The VACUETTE® Urine Collection System enables a closed, sterile urine collection that in turn will maintain the sample quality and guarantee the efficiency of the testing procedure. The product line includes tubes, cups, sets and various accessories.

URINE COLLECTION SYSTEM

/	VACUETTE® Urine Tubes	154
	VACUETTE® Z Urine No Additive Tubes	155
	VACUETTE® Urine CCM Tubes	156
/	Urine Sets	157
	Urine Sets	158
	Urine CCM Sets	158
/	Accessories Urine Collection	160
	Urine Transfer Devices	161
	Urine Cups	161
	24h Urine Containers	162



/ You can find more information about this product in our web-shop.



VACUETTE® URINE TUBES

VACUETTE® Urine Tubes are used for collection, transportation and/or analysis. They are sterile, leak-proof and virtually unbreakable.

VACUETTE® Urine Tubes are available with either a round or conical base. Round-based tubes are used for chemical urine tests, while tubes with a conical base are used for microscopic examinations of urine sedimentation.

Our tubes are also available with

or without preservatives. In the case of the tubes with preservative, VACUETTE® Urine CCM tubes help stabilize the sample, which is important when this one cannot be analyzed within 2 hours of collection. The sample collected in it should be inverted 6-8 times following collection in order to achieve a homogenous mixing.

All VACUETTE® Urine Tubes are also available with a pre-barcoded label.

- / Round base or conical base
- / Without additive
- / With preservative (CCM)
- / Available pre-barcoded



IVD


STERILE

VACUETTE®

Z Urine No Additive Tubes

Z Urine No Additive Tubes are used as collection and/or transport tubes in chemical urine analysis, urine test strips, for microbiological studies and urine microscopy. It is recommended to perform the analysis within a short time after collection.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Bottom shape	Label	Barcode	Qty. inner / outer
454141	3 ml	yellow	yellow	PREMIUM	13 x 75	round	paper	no	50 / 1,200
484533	4 ml	yellow	yellow	PREMIUM	13 x 75	round	paper	yes	50 / 1,200
456007	6 ml	yellow	yellow	PREMIUM	13 x 100	round	paper	no	50 / 1,200
476081	6 ml	yellow	yellow	PREMIUM	13 x 100	round	transparent	no	50 / 1,200
486027	6 ml	yellow	yellow	PREMIUM	13 x 100	round	paper	yes	50 / 1,200
456065	6 ml	yellow	yellow	non-ridged	13 x 100	round	paper	no	50 / 1,200
485523	7 ml	yellow	yellow	non-ridged	16 x 100	round	paper	yes	50 / 1,200
455227	7 ml	yellow	yellow	non-ridged	16 x 100	round	paper	no	50 / 1,200
455225	8 ml	yellow	black	non-ridged	16 x 100	round	paper	no	50 / 1,200
455028	9 ml	yellow	yellow	non-ridged	16 x 100	conical	paper	no	50 / 1,200
485507	9 ml	yellow	yellow	non-ridged	16 x 100	conical	paper	yes	50 / 1,200
495028	9 ml	yellow	yellow	non-ridged	16 x 100	conical	BiBarcode (paper)	yes	50 / 1,200
485510	10 ml	yellow	yellow	non-ridged	16 x 100	round	paper	yes	50 / 1,200
455007	10 ml	yellow	yellow	non-ridged	16 x 100	round	paper	no	50 / 1,200
495007	10 ml	yellow	yellow	non-ridged	16 x 100	round	BiBarcode (paper)	yes	50 / 1,200

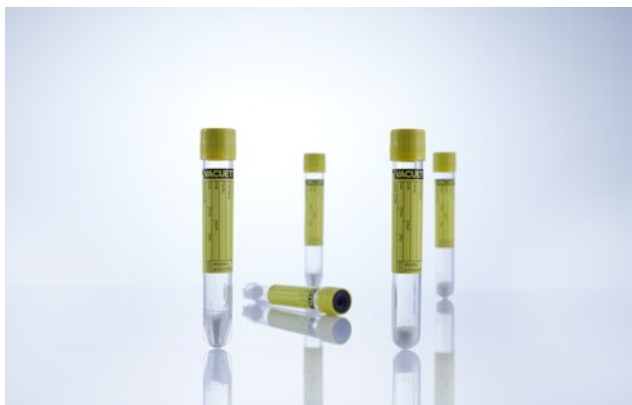

IVD

STERILE

VACUETTE®

Urine CCM Tubes

The VACUETTE® Urine CCM Tubes contain a powdered stabilizer that keeps the bacteria and yeasts in the sample stable at room temperature (20 - 25°C) for up to 48 hours. The CCM tube can thus win time for carrying out the analysis.

Item no.	Nominal volume	Cap color	Ring color	Thread type	Tube size	Bottom shape	Label	Barcode	Qty. inner / outer
454486	4 ml	● yellow	● black	non-ridged	13 x 75	round	paper	no	50 / 1,200
484532	4 ml	● yellow	● black	non-ridged	13 x 75	round	paper	yes	50 / 1,200
456284	6 ml	● yellow	● black	non-ridged	13 x 100	round	paper	no	50 / 1,200
496510	6 ml	● yellow	● black	non-ridged	13 x 100	round	BiBarcode (paper)	yes	50 / 1,200
455243	9 ml	● yellow	● black	non-ridged	16 x 100	conical	paper	no	50 / 1,200
455052	10 ml	● yellow	● black	non-ridged	16 x 100	round	paper	no	50 / 1,200
485511	10 ml	● yellow	● black	non-ridged	16 x 100	round	paper	yes	50 / 1,200



/ You can find more information about this product in our web-shop.



URINE SETS

Urine tubes and the Urine Transfer Device from Greiner Bio-One are also available in a practical set.

VACUETTE® Urine Tubes are available in a set with the Urine Transfer Device to provide more convenience for daily routine work. With the urine sets from Greiner Bio-One, you have all the products you need at hand at all times.

The set is available with round bottom tubes or conical bottom tubes. You can also choose between a set with tubes without additives or a set with tubes with preservative (CCM).

The Urine Transfer Device included in the set allows you to transfer the sample into the tubes easily and safely.

- / VACUETTE® Urine Tubes or VACUETTE® Urine CCM Tubes
- / Round base or conical base
- / Urine Transfer Device



IVD

Urine Sets

The VACUETTE® Urine Tube is also available in a practical set: VACUETTE® Urine Tube + Urine Transfer Device.

Item no.	Description	Nominal volume	Bottom shape	Qty. inner / outer
453070	Urine Set 16x100 (VACUETTE® Urine Tube Item no. 455028 + Urine Transfer Device Item no. 450751) single packed	9 ml	conical	- / 300



IVD

Urine CCM Sets

The VACUETTE® CCM Urine Tube is also available in a practical set: VACUETTE® CCM Urine Tube + Urine Transfer Device. The set is available with both the round bottom tube and the conical bottom tube.

Item no.	Description	Nominal volume	Bottom shape	Qty. inner / outer
453031	Urine CCM Set 13x75 (VACUETTE® Urine CCM Tube Item no. 454486 + Urine Transfer Device Item no. 450751) single-packed	4 ml	round	1 / 300
453032	Urine CCM Set 16x100 (VACUETTE® Urine CCM Tube Item no. 455243 + Urine Transfer Device Item no. 450751) single-packed	9 ml	conical	- / 300
453034	Urine CCM Set 16x100 (VACUETTE® Urine CCM Tube Item no. 455052 + Urine Transfer Device Item no. 450751) single-packed	10 ml	round	- / 300

Item no.	Description	Nominal volume	Bottom shape	Qty. inner / outer
483070	Urine CCM Set 16x100 (VACUETTE® Urine CCM Tube Item. no. 485511 + Urine Transfer Device Item.no. 450751) with barcode label	10 ml	round	- /300



/ You can find more information about this product in our web-shop.




ACCESSORIES URINE COLLECTION

The Greiner Bio-One Urine Line is rounded off by a comprehensive range of accessories.

In addition to the VACUETTE® Urine Tubes and sets, Greiner Bio-One also includes urine cups in its range, which are available with integrated transfer device

to ensure clean sample transfer. The range of accessories also includes 24 Hours Urine Containers with or without integrated transfer device as well as the Urine Transfer Device itself in the long or short version – according to your needs.

- / Wide range of accessory products
- / Urine cups, sterile or non-sterile
- / 24 Hours Urine Container
- / Products available with or without integrated transfer device
- / Urine Transfer Device

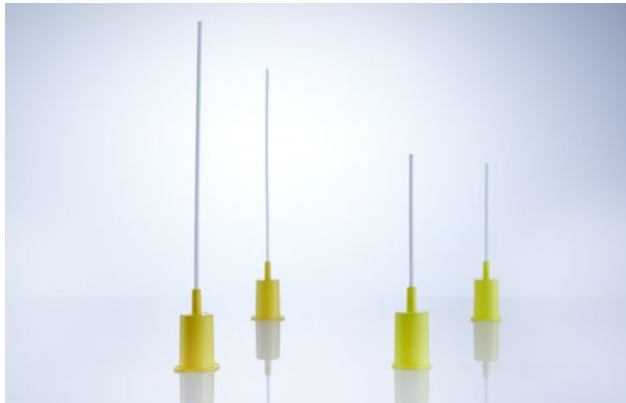

IVD

Urine Transfer Devices

The Urine Transfer Device allows easy and hygienic transfer of sample material from a collection transfer container. The product is available in two lengths.

Item no.	Description	Qty. inner / outer
450251	short	50 / 600
450252	long	50 / 600
450253V1	short, single-packed, sterile	- /400
450264	short, single-packed, non-sterile	- /400



IVD

Urine Cups

Urine cups from Greiner Bio-One come with a transfer device already integrated in the cap. This features gives you the advantage of a closed system, ensuring a clean and hygienic sample transfer into a tube.

Cap design: screw cap with integrated transfer device

Item no.	Description	Cap color	Qty. inner / outer
724321	100 ml, bulk-packed, sterile	● yellow	- /300
724322	100 ml, bulk-packed, non-sterile	● yellow	- /300
724325	100 ml, 10 pcs./bag, sterile	● yellow	10/250



24h Urine Containers

The 24h urine collection container is available as a 2.7 liter version or as VACUETTE® URI-PLUS with a 3 liter volume. VACUETTE® URI-PLUS has an integrated transfer device for even more hygienic sample collection.

Graduation: yes

Item no.	Description	Product color	Cap color	Cap design	Qty. inner / outer
724315	24h urine container	● brown	○ white	screw cap	- /34
724316	VACUETTE® URI-PLUS	● brown	● blue	screw cap with integrated transfer device	- /36



Greiner Bio-One's product range for molecular biology and virology covers tubes for storage of the sample and transport to its laboratory destination.

MOLECULAR BIOLOGY

/	VACUETTE® Virus Stabilization Tube ...	166
	VACUETTE® Virus Stabilization Tube	167



/ You can find more information about this product in our web-shop.



VACUETTE®

VIRUS STABILIZATION TUBE

In the VACUETTE® Virus Stabilization Tube from Greiner Bio-One, viruses are safely contained and can be used for the detection of SARS CoV-2, Influenza A and Influenza B.

With the VACUETTE® Virus Stabilization Tube, Greiner Bio-One has been able to make a significant contribution to virus detection since the outbreak of the

COVID 19 pandemic in 2020. Respiratory viruses remain a major health concern worldwide, and ensuring proper pre-analytical processes is critical to achieving accurate and dependable diagnostic outcomes. Greiner Bio-One offers products that support safe transport (see chapter «Sample Transport»)

/ For the detection of SARS CoV-2, Influenza A and Influenza B



VACUETTE® Virus Stabilization Tube

VACUETTE® Virus Stabilization Tubes are PET tubes without vacuum and sterile inside. They contain a phosphate-buffered saline solution in various volumes. The tube is intended for storage and transport of viral sample material.

Item no.	Nominal volume	Cap color	Thread type	Tube size	Barcode	Qty. inner / outer
456162	2 ml	● red	PREMIUM	13 x 100	no	50 / 1,200
456161	3 ml	● red	PREMIUM	13 x 100	no	50 / 1,200

NUMERICAL INDEX

370267B.....	147	450252.....	161	454081.....	49	455040.....	54	456284.....	156
371516.....	145	450253V1.....	161	454083.....	50	455045.....	53	456287A.....	50
371517.....	145	450263.....	115	454083A.....	50	455051.....	51	456305A.....	50
371518.....	145	450264.....	161	454086.....	55	455052.....	156	456309.....	47
371519.....	145	450330V1.....	133	454087.....	54	455057.....	45	456310A.....	47
371522.....	145	450331V1.....	133	454088.....	57	455057A.....	45	456313.....	47
371523.....	145	450332.....	133	454089.....	49	455071.....	45	458711.....	126
371527.....	145	450333.....	133	454092.....	44	455071A.....	45	458712.....	126
450040.....	109	450334.....	133	454096.....	44	455072.....	46	458713.....	126
450041.....	109	450335.....	133	454098.....	44	455072A.....	46	458714.....	126
450042.....	109	450338.....	134	454141.....	155	455083.....	50	458715.....	126
450043.....	109	450339.....	134	454204.....	44	455083A.....	50	459000.....	119
450064.....	111	450343.....	128	454205A.....	45	455084.....	49	459001.....	119
450065.....	111	450373.....	111	454209.....	53	455092.....	44	459010.....	119
450066.....	111	450417.....	147	454214.....	45	455099A.....	46	459011.....	119
450070.....	121	450418.....	147	454214A.....	45	455222.....	49	459020.....	119
450071.....	109	450427.....	99	454217.....	55	455225.....	155	459021.....	119
450072.....	109	450428.....	99	454221.....	59	455227.....	155	472001.....	139
450073.....	109	450429.....	99	454222.....	55	455229.....	50	472010.....	139
450075.....	109	450438.....	99	454225A.....	45	455229A.....	50	472011.....	139
450076.....	109	450439.....	99	454228.....	45	455243.....	156	472015.....	139
450077.....	109	450521.....	98	454228A.....	45	456001.....	57	472016.....	139
450081V1.....	95	450522.....	98	454235.....	54	456004A.....	46	472017.....	139
450082V1.....	95	450523.....	98	454236.....	44	456005A.....	46	472020.....	141
450083V1.....	95	450524.....	98	454237.....	49	456006A.....	46	472022.....	141
450084V1.....	95	450530.....	79	454241.....	57	456007.....	155	472023.....	141
450085V1.....	96	450531.....	79	454243.....	45	456010.....	45	472024.....	141
450086V1.....	96	450532.....	79	454243A.....	45	456010A.....	45	472030.....	141
450089.....	111	450533.....	75	454246.....	53	456018.....	46	472040.....	139
450090V1.....	95	450534.....	75	454247A.....	51	456018A.....	46	474036.....	55
450091V1.....	95	450535.....	77	454253A.....	45	456019A.....	46	474064.....	42
450092V1.....	95	450536.....	77	454260A.....	45	456023.....	53	474065.....	42
450095V1.....	95	450537.....	77	454263A.....	45	456036.....	55	474071A.....	46
450096V1.....	95	450539.....	73	454287.....	55	456038.....	55	474083A.....	50
450099V1.....	95	450540.....	81	454297.....	59	456043.....	53	474084.....	49
450121.....	92	450541.....	81	454319.....	45	456048.....	44	474086.....	55
450122.....	93	450545.....	79	454319A.....	45	456052.....	64	474087.....	54
450123.....	93	450546.....	79	454321.....	41	456058.....	54	474092.....	44
450125.....	93	450547.....	79	454322.....	41	456065.....	155	474214.....	46
450126.....	93	450548.....	75	454325.....	41	456069.....	45	474214A.....	46
450127.....	92	450549.....	75	454327.....	41	456069A.....	45	474321.....	41
450128.....	92	450550.....	77	454332.....	41	456071.....	46	474327.....	41
450129.....	92	450551.....	77	454334.....	41	456071A.....	46	476010.....	46
450130.....	92	450552.....	81	454385.....	42	456073.....	46	476010A.....	46
450132.....	92	450568.....	75	454410.....	53	456073A.....	46	476071.....	46
450159.....	111	453031.....	158	454411.....	53	456080.....	62	476071A.....	46
450160V1.....	96	453032.....	158	454420A.....	46	456083.....	50	476081.....	155
450161V1.....	96	453034.....	158	454483A.....	50	456083A.....	50	476083.....	50
450181.....	117	453070.....	158	454486.....	156	456084.....	49	476092.....	44
450182.....	96	454001.....	57	454489.....	50	456085.....	57	480539.....	73
450183.....	96	454008.....	50	454489A.....	50	456087.....	50	480546.....	79
450184.....	96	454008A.....	50	454491A.....	50	456087A.....	50	480547.....	79
450185.....	96	454020.....	53	454493A.....	50	456088.....	49	480548.....	75
450191V1.....	95	454021.....	55	454510.....	60	456089.....	44	480549.....	81
450209.....	115	454023.....	53	454511.....	60	456092.....	44	480550.....	77
450210.....	116	454024.....	53	454513.....	60	456093.....	64	480560.....	77
450216.....	115	454028.....	45	454514.....	60	456097.....	64	483070.....	159
450225.....	116	454028A.....	45	454591.....	47	456161.....	167	484030A.....	51
450226V1.....	90	454036.....	55	454591A.....	47	456162.....	167	484031A.....	51
450228V1.....	90	454046A.....	50	454592.....	47	456211.....	54	484507.....	53
450230.....	89	454047.....	53	454593.....	47	456214A.....	46	484508.....	41
450235V1.....	90	454067.....	45	455001.....	57	456216.....	50	484509.....	54
450237V1.....	89	454067A.....	45	455007.....	155	456216A.....	50	484510.....	55
450239V1.....	90	454071.....	45	455021.....	54	456234.....	46	484511.....	55
450241.....	115	454071A.....	45	455028.....	155	456234A.....	46	484512.....	54
450251.....	161	454078A.....	45	455036.....	55	456275.....	62	484513.....	46

484513A	46	VP164511	104
484514	46	VP174501	104
484514A	46	VP174511	104
484515	44	VP183201	104
484517	41	VP183211	104
484518	41	VP184501	104
484520	59	VP184511	104
484522	49	VP202501	104
484523	50	VP202511	104
484523A	50	VP203201	104
484525	44	VP203211	104
484526	53	VP222501	104
484527	59	VP222511	104
484528	59	VP241901	104
484530	57	VP241911	104
484531	51	VW144501	103
484532	156	VW144511	103
484533	155	VW164501	104
484586	53	VW164511	103
484589	49	VW174501	104
485222	49	VW174511	103
485501	55	VW183201	104
485502	46	VW183211	103
485506	57	VW184501	104
485507	155	VW184511	103
485510	155	VW203201	104
485511	156	VW203211	103
485513	50	VW222501	104
485514	46	VW222511	103
485514A	46	VW241911	103
485523	155		
486004A	46		
486005A	46		
486006A	46		
486007A	46		
486027	155		
486500	46		
486502	55		
486503	46		
486504	44		
486505	41		
486507	50		
486508	49		
486509	47		
495007	155		
495028	155		
496510	156		
724315	162		
724316	162		
724321	161		
724322	161		
724325	161		
800101	139		
800102	139		
800105	139		
800110	139		
840050	123		
840051	123		
840052	124		
840053	123		
890001	151		
890004	151		
NW241901	103		
NW261901	103		
VP144501	104		
VP144511	104		
VP164501	104		

making a difference

making a difference

www.gbo.com

GREINER BIO-ONE GMBH
KREMSMÜNSTER, ÖSTERREICH

PHONE +43 7583 6791-0
FAX +43 7583 6318
E-MAIL office@at.gbo.com



GREINER BIO-ONE IS A GLOBAL PLAYER.
FIND THE CONTACT DETAILS OF YOUR
LOCAL PARTNER ON OUR WEBSITE.



This product information is addressed exclusively to healthcare professionals. Devices of Greiner Bio-One are to be used by properly trained healthcare professionals only in accordance with the relevant Instructions for Use (IFU). For a listing of indications, contraindications, precautions and warnings, please refer to the Instructions for Use which accompanies each product or is available for download from our website at www.gbo.com (Download Center). For more information contact your local Greiner Bio-One sales representative or visit our website.

All information is provided without guarantee despite careful processing. Any liability, warranty or guarantee of Greiner Bio-One GmbH is excluded. All rights, errors and changes are reserved. If not stated otherwise, Greiner Bio-One GmbH has all copyrights and/or other (user-)rights in this documents, in particular to signs such as the mentioned (word-picture-) brands and logos. Any use, duplication or any other use of the rights of Greiner Bio-One GmbH is expressly prohibited.
Media owner: Greiner Bio-One GmbH / **Manufacturer:** [Samson Druck GmbH / 5581 St. Margarethen]

980042 en [rev.14 02.2026] - Not valid for the US market


greiner
BIO-ONE